


# Ultrafast Fiber Lasers

The Next Generation – Reliable, Robust, Flexible



*A Passion for Precision.*

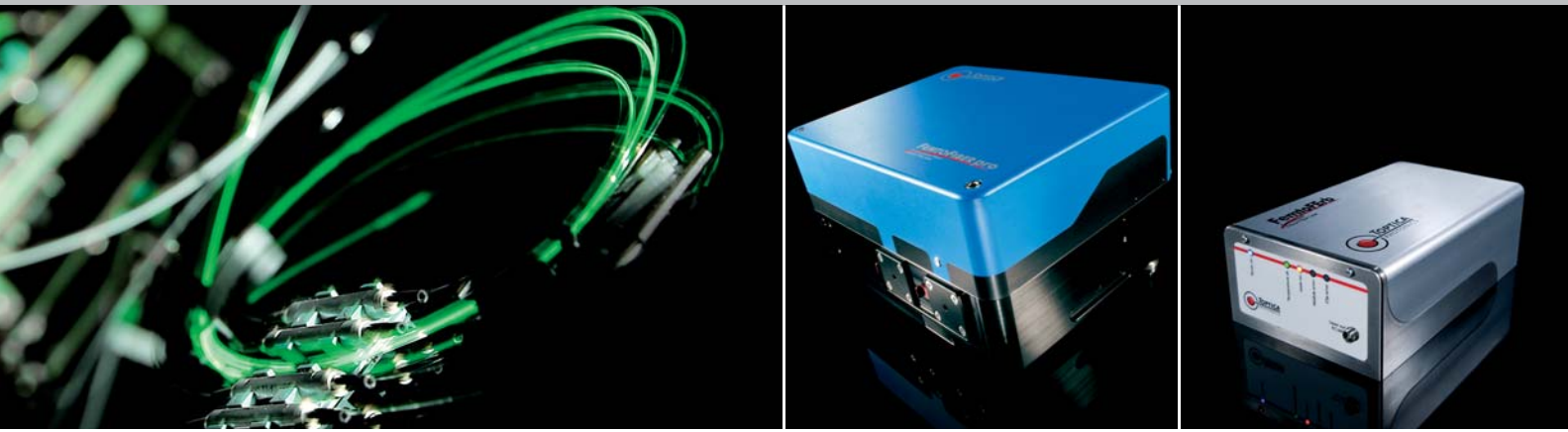



# FemtoFiber® Lasers — The Next Generation of Ultrafast Technology



## Contents

Ultrafast technology for professional applications .....	4 - 5
Fiber laser technology .....	6 - 7
FemtoFiber pro — modularity .....	8 - 11
IR systems (1560 nm) .....	12
NIR systems (780 nm) .....	13
SCIR systems (980 nm - 2200 nm) .....	14
UCP systems (980 nm - 1400 nm, < 25 fs) .....	15
TVIS systems (488 nm - 640 nm) .....	16
TNIR systems (830 nm - 1100 nm) .....	17
Options (VAR, LRC, M40, AMP) .....	18
FemtoFiber pro — control unit .....	19
Industrial fiber lasers .....	
iChrome (488 nm - 640 nm) .....	20
PicoFYb (1030 / 1064 nm), FemtoFErb (1560 nm) ...	21

FemtoFiber pro

iChrome

FemtoFErb / PicoFYb



Ultrafast technology has seen an unpredictable success ever since it was introduced. Today physicists are still pushing the limits, e.g. the pulse duration down to a single optical cycle and the peak power up to tera- or even petawatt. On the other hand, many promising applications have emerged outside the laser laboratories. Mainly the high peak power and ultrashort pulse duration make ultrafast lasers very attractive. Both give rise to nonlinear effects and open new paths in engineering and scientific research.

Independent of the application, the key for successful integration of ultrafast lasers is the reduction of complexity. Many years ago, only laser experts could handle such sensitive devices. Today, one expects an ultrafast laser to work by pushing a single button. TOPTICA's answer to these expectations resulted in several matured products: The FemtoFiber pro series, iChrome, PicoFYb and the FemtoFErb.

First of all, these fiber laser systems are based on polarization-maintaining fibers. Reliable mode-locking is achieved with a saturable absorber mirror (SAM). Only high quality components with extraordinary long lifetimes are built in. These comply with telcordia specifications and are suitable for 24/7 operation. The heart of each laser – namely the master oscillator and the amplifier – is completely fiber-based which makes alignment redundant and the systems extremely robust.


These and other advantages helped TOPTICA's ultrafast fiber lasers to become the source of choice for the following applications: Life sciences, time domain terahertz spectroscopy, attoscience and material processing among others.



## Ultrafast Technology for Professional Applications

### Life sciences

Modern biomedical research often involves advanced laser microscopy. Tunable lasers account for the variety of fluorescent labels. Pulsed lasers enable fluorescent lifetime imaging microscopy (FLIM). Ultrafast lasers open the door for noninvasive methods relying on multiphoton excitation such as SHG, THG or CARS microscopy. While the techniques become more advanced, the

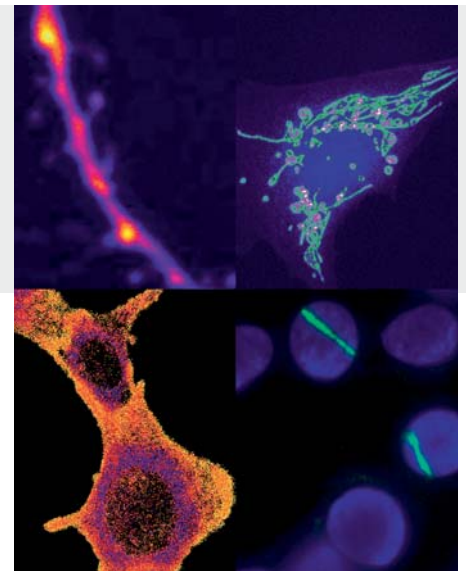


employed lasers should become easier to handle. That's why TOPTICA's fiber lasers are the product of choice.

Reliable and comfortable in handling, they allow the users to focus on their scientific task. Especially the iChrome is a fully automated and fiber coupled laser designed for confocal microscopy in the visible range (488 nm - 640 nm). On the other hand, the unique flexibility makes the FemtoFiber pro a versatile tool: Not only picosecond laser pulses that are tunable from 488 nm to 640 nm but also femtosecond pulses at 780 nm or tunable from 980 nm to 1400 nm can be combined in a single laser system.

## 2-photon polymerization

Miniaturized devices will change our life: Maybe someday micro-meter sized robots will travel through our blood circuit system and repair our cells. Today, nanotubes, nanowires, etc. are fabricated in nanotechnology centers around the world. For example, photonic crystals are produced by two-photon polymerization (2PP): The exciting property of such

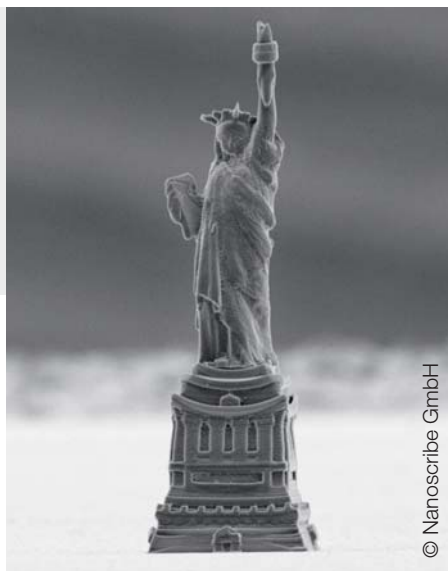

crystals is their photonic band gap which prohibits that light of certain wavelengths or polarizations can propagate through such a device. This band gap can be designed by the type and dimensions of the structure itself which enables the designer to directly access optical properties.

TOPTICA delivers the fiber laser with highest peak power at 780 nm that can be easily integrated in a commercial 2PP microscope. This enables the user to write structures with dimensions as big as 300  $\mu\text{m}$  along the optical axis — the highest reported with fiber lasers.

## Time-domain terahertz

Terahertz technology features the chemical identification of substances through opaque materials such as cloth and synthetics, e.g for quality control and security tasks. Broadband terahertz radiation is generated with femtosecond laser pulses: Either they are focused on antenna devices such as photoconductive switches based on GaAs or InP, or they undergo optical rectification




in nonlinear crystals, such as GaP, ZnTe or DAST. For high power and large bandwidth of the THz radiation peak power and pulse width of the laser are crucial.

TOPTICA offers two ideal laser sources for THz generation: The FemtoFiber pro and the FemtoFERb. Both systems provide outstanding peak power levels at pulse widths well below 100 fs. The FemtoFiber pro NIR is the ideal solution for research labs, thanks to its beneficial switching possibility between 780 and 1560 nm. Moreover, fast optical sampling of the THz signal can be accomplished without moving parts employing the ECOPS option. The FemtoFERb as an ultra-compact fiber-coupled system on the other hand is designed for easiest OEM integration.

## Material processing

The trend in laser micro-machining goes from nanosecond to picosecond laser pulses. Cold ablation without thermal load is the key for unprecedented precision. Furthermore, fiber lasers will replace solid state lasers in this market. Hands-off and maintenance-free operation, robustness and compactness, low running costs and better thermal properties are but only a few

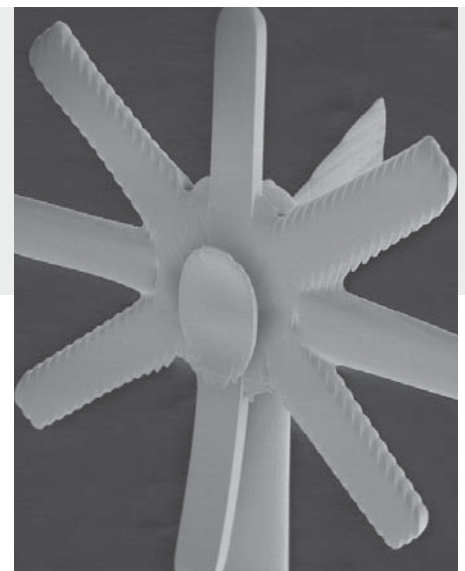

of the advantages over free-space, bulky laser resonators.

TOPTICA's PicoFYb is specially designed for OEM integration in laser amplifiers. Most stable laser operation is guaranteed in respect to amplitude noise, pulse duration, frequency jitter and reliable mode-locking. The compact laser head includes all driver electronics and delivers the laser pulses in a polarization-maintaining fiber. The typical wavelengths are 1030 nm and 1064 nm.



# Fiber Laser Technology

## FemtoFiber pro

### Mode of operation


The engine of all TOPTICA's fiber lasers is a SAM mode-locked ring oscillator. This configuration – protected by an optical isolator – is known as most robust against back reflections of any kind. Running in the solitonic regime, the oscillator operates in a well defined state. Special TOPTICA design ensures full suppression of parasitic pedestals (Kelly-bands). The soliton pulses are carefully amplified to very high peak powers in a subsequent core-pumped

fiber amplifier. This MOPA (master oscillator power amplifier) platform is all-fiber and based completely on polarization-maintaining fibers. In contrast to single-mode fiber compressors, the motorized silicon prism compressor preserves the smooth spectral shape. This results in nearly pedestal-free sub-100 fs pulse shapes in the time domain. Optional modules convert the center wavelength into supercontinua, tunable visible or NIR laser light.




### PM fiber technology

Polarization-maintaining (PM) fibers preserve the linear state of polarization under all circumstances: Temperature, mechanical stress, humidity or other environmental influences can't disturb the stable laser operation.

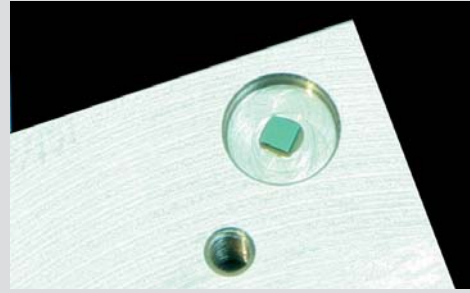

### SAM mode-locking technology

The saturable absorber mirror (SAM) ensures self-starting and reliable cw-mode-locking. This passive device efficiently suppresses Q-switching operation. Only intense solitonic laser pulses are stabilized and supported by the SAM.

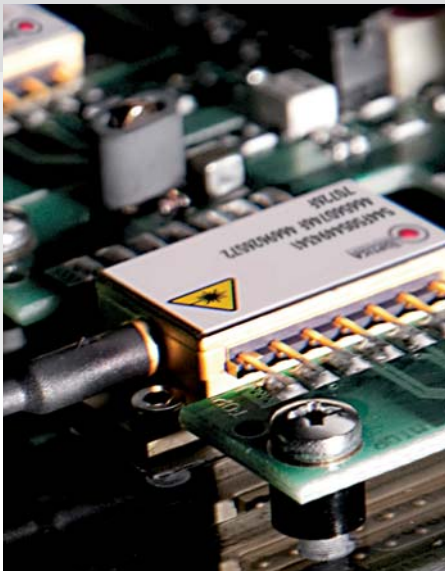

### Robust 24/7 operation


- PM fiber assembly
- All-fiber: Alignment free
- Telcordia proved components

### Turnkey

- Reliable SAM technology
- Self-starting mode-locking

### Price and system size

- Low cost of ownership & running costs
- Compact: Letter / A4 format footprint
- Neither air nor water cooling

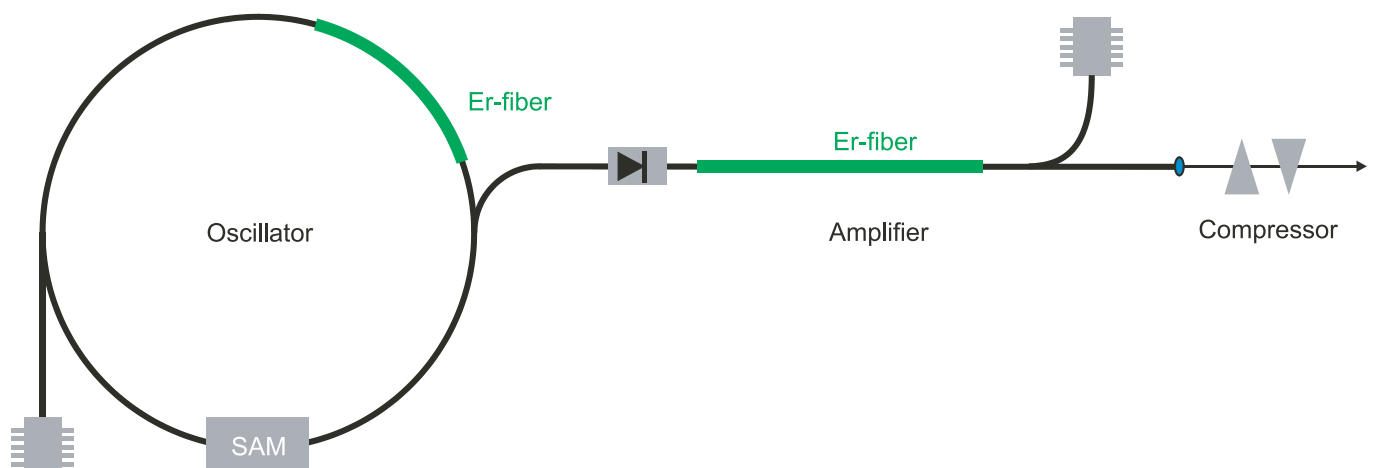



## Modularity

### FemtoFiber pro




#### Versatile laser system

An all-fiber master oscillator power amplifier (MOPA) concept is the basis for TOPTICA's ultrafast fiber lasers. Either Ytterbium or Erbium doped gain fibers support ultrashort laser pulses at 1030 nm, 1064 nm or 1560 nm. Depending on the application, femtosecond or picoseconds pulse widths are custom tailored. Optional modules extend the available wavelengths:


Second harmonic generation adds 780 nm. Highly nonlinear fibers generate an octave spanning supercontinuum (980 nm - 2200 nm). Subsequent second harmonic generation results in tunable laser light in the visible range (488 nm - 640 nm). Fiber delivery makes integration easy. Finally, customization completes our idea of flexibility.



# Outstanding Performance

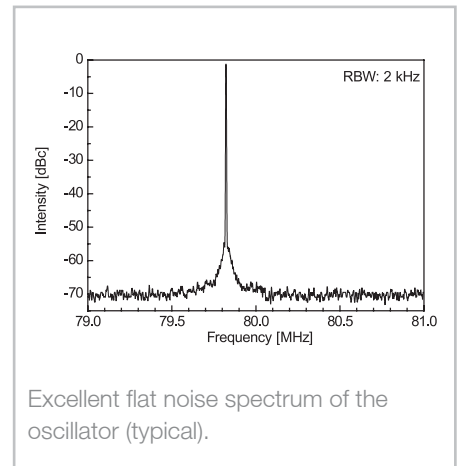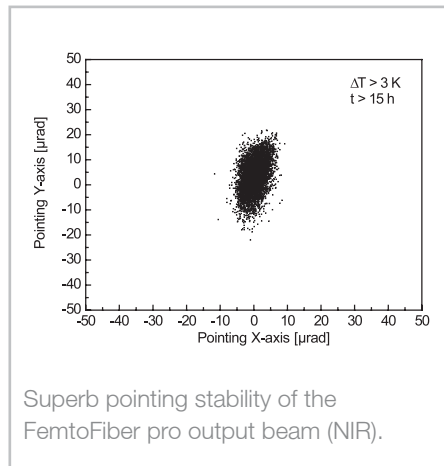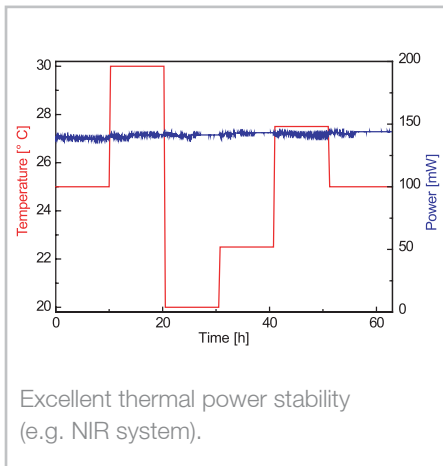

# Challenging Applications

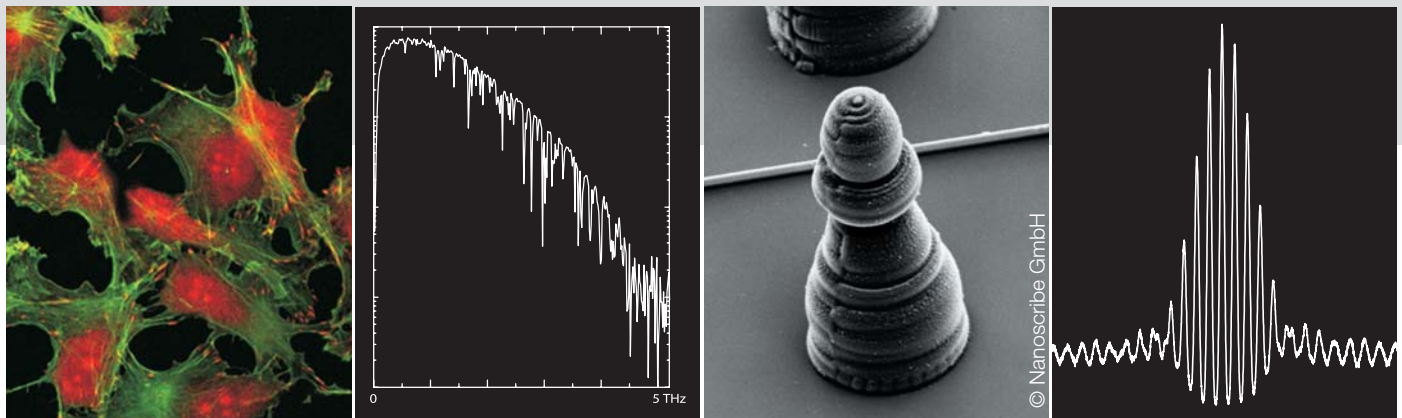

## Life sciences

### Laser requirements:

- Tunability
- Ultrashort pulses
- Multi color

### Recommended products:

- iChrome
- FemtoFiber pro TVIS
- FemtoFiber pro TNIR
- FemtoFiber pro NIR

## Time-domain terahertz

### Laser requirements:

- < 100 fs pulse width
- 1560 nm or 780 nm

### Recommended products:

- FemtoFiber pro IR
- FemtoFiber pro NIR
- FemtoFErb 1560

## 2-photon polymerisation

### Laser requirements:

- High laser peak power
- Short pulses

### Recommended products:

- FemtoFiber pro NIR






## Attoscience

### Laser requirements:

- CEO-free laser pulses
- Mid-IR radiation for parametric amplification

### Recommended products:

- FemtoFiber pro IR + FemtoFiber pro SCIR
- FemtoFiber pro UCP

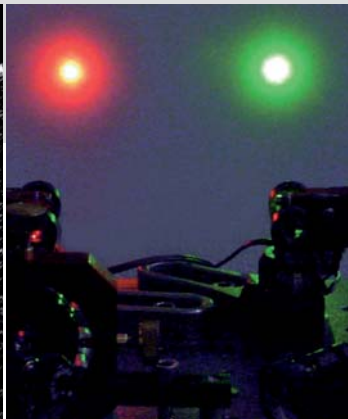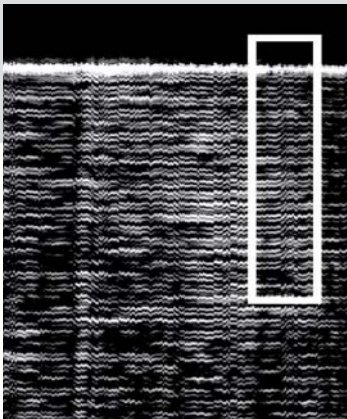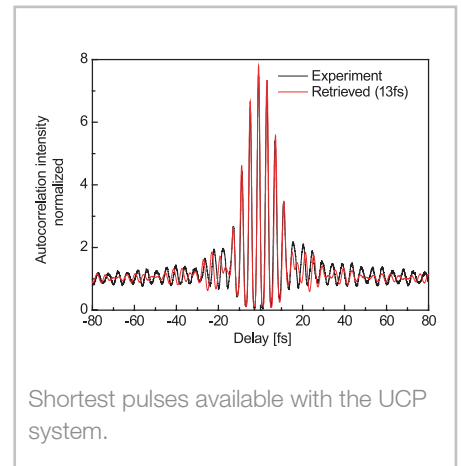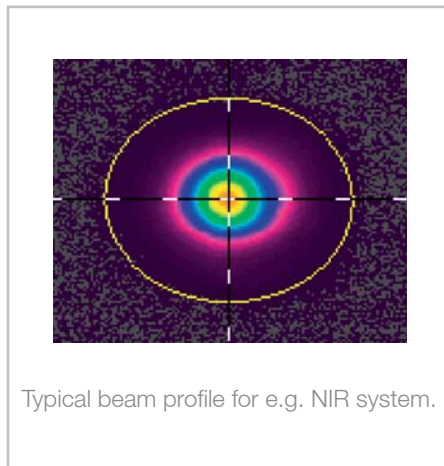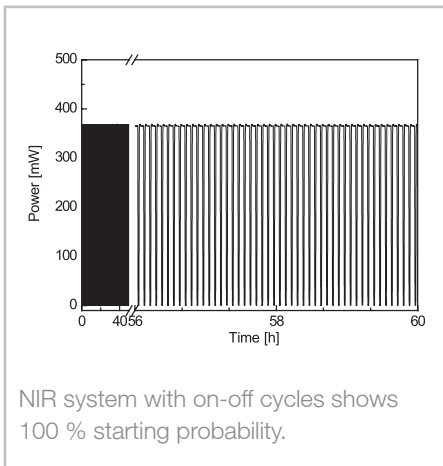

### Optical coherence tomography

- Laser requirements:**
- Very broad spectrum
  - Supercontinuum

- Recommended products:**
- FemtoFiber pro IR
  - FemtoFiber pro SCIR

### Pump-probe spectroscopy

- Laser requirements:**
- Phase coherent dual beam system
  - Independently configured outputs

- Recommended products:**
- FemtoFiber pro with multi beam option

### Material processing

- Laser requirements:**
- Picosecond pulses
  - 1030 nm or 1064 nm

- Recommended products (seeders):**
- PicoFYb 1030 nm AMP
  - PicoFYb 1064 nm AMP

### ECOPS

- Laser requirements:**
- Synchronization of the repetition rate
  - Phase-locked loop electronics

- Recommended products:**
- FemtoFiber pro with VAR option
  - FemtoFiber pro laser repetition rate control (LRC)

## FemtoFiber pro specifications

Models	Wavelength	Pulse width	Average power	Page
FemtoFiber pro <b>IR</b>	1560 nm	< 100 fs	> 350 mW	12
FemtoFiber pro <b>NIR</b>	780 nm	< 100 fs	> 140 mW	13
FemtoFiber pro <b>SCIR</b>	980 - 2200 nm	< 100 fs	> 150 mW	14
FemtoFiber pro <b>UCP</b>	980 - 1400 nm	< 25 fs	> 30 mW	15
FemtoFiber pro <b>TVIS</b>	488 - 640 nm	< 1 ps	1 – 10 mW	16
FemtoFiber pro <b>TNIR</b>	830 - 1100 nm	< 200 fs	1 – 5 mW	17


## Options

VAR	Modulation of the pulse repetition frequency Piezo-transducer resonance frequency > 5 kHz Repetition frequency tuning range > 200 kHz	18
LRC	Phase-locked loop electronics RMS jitter t.b.d. Electronics rack with power supply USB interface and control software	18
M40	Oscillator design with 40 MHz repetition rate (some specifications may change)	18
AMP	System without oscillator for multi beam systems FC/APC fiber input	18


## General specifications

Repetition rate	80 MHz standard
Output coupling	Free space
Beam shape	TEM <sub>00</sub>
Polarization	Linear, horizontal
Laser head dimensions (w x d x h)	<b>IR, NIR, SCIR:</b> 280 x 229 x 151 mm <sup>3</sup> <b>UCP, TVIS, TNIR:</b> 280 x 458 x 151 mm <sup>3</sup>
Control unit dimensions (w x h x d)	235 x 315 x 140 mm <sup>3</sup>
Line input	90 – 260 VAC, 47 – 63 Hz
PC interface	Ethernet, USB, RS 232


Specifications are subject to change without further notice

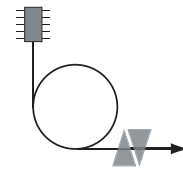


FemtoFiber pro **IR** — System unit consisting of mode-locked Erbium fiber oscillator with integrated fiber amplifier and dispersion control module.

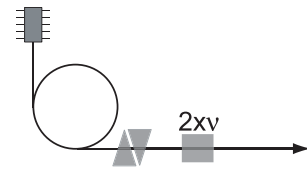


FemtoFiber pro **NIR** — System unit with additional second harmonic generation.

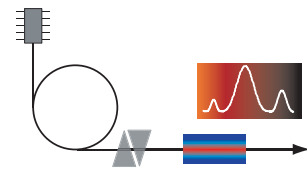


FemtoFiber pro **SCIR** — System unit with additional highly nonlinear fiber generating octave-spanning supercontinuum in the infrared.

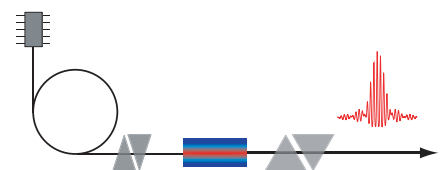


FemtoFiber pro **UCP** — Laser system generating < 25 fs pulses, continuously tunable from 980 – 1400 nm.

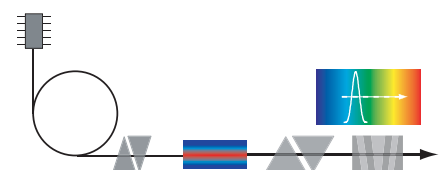

FemtoFiber pro **TVIS** — Laser system with ultrashort pulsed beam, continuously tunable from 488 – 640 nm.

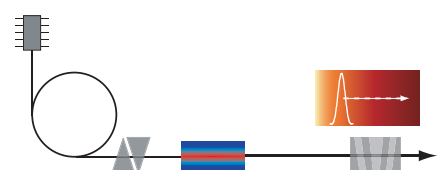

FemtoFiber pro **TNIR** — Laser system with ultrashort pulsed beam, continuously tunable from 830 – 1100 nm.

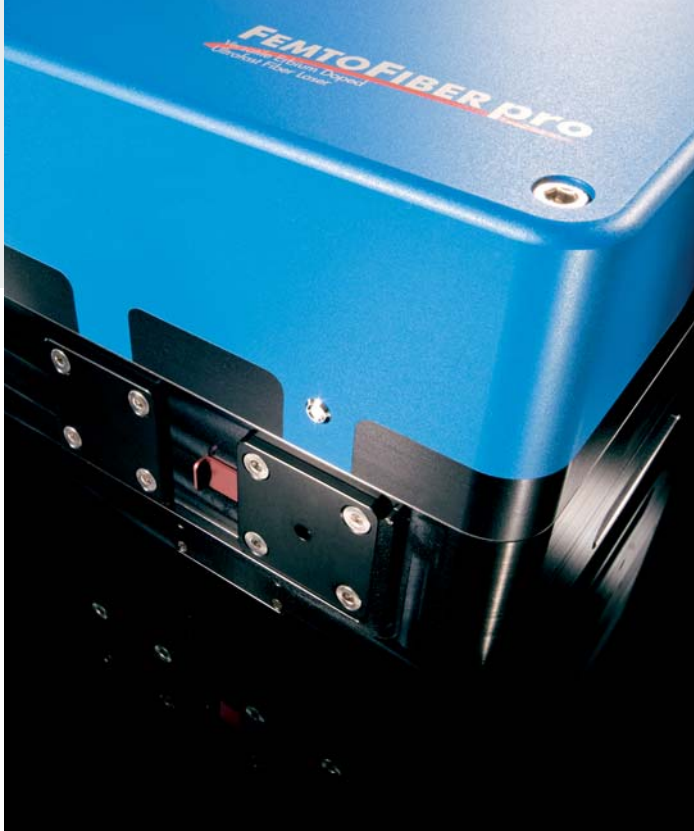

# FemtoFiber pro

**IR**  
InfraRed

## Key features




- Fiber laser system with SAM mode-locked oscillator and core-pumped high power amplifier in a single box
- Control unit with laser drivers, power PC control, communication interfaces and power supply
- Motorized prism compressor for pulse optimization
- Free-beam output
- Secondary oscillator fiber output (FC/APC connector)
- High bandwidth monitor output (SMA connector)

## Key specifications

Center wavelength	<b>1560 nm</b>
Laser output power	> 350 mW
Pulse width	< 100 fs
Repetition rate	80 MHz standard
Linear polarization	> 95 % (horizontal)
Beam shape	TEM <sub>00</sub> , M <sup>2</sup> < 1.2
Beam size (1/e <sup>2</sup> )	Typ. Ø 3.5 mm
Beam divergence	< 2 mrad
Output coupling	Free space
Laser head dimensions	280 x 229 x 151 mm <sup>3</sup> (w x d x h)
Laser head weight	< 10 kg
Control unit dimensions	235 x 315 x 140 mm <sup>3</sup> (w x d x h)
Control unit weight	< 4.5 kg
Line input	90 – 260 VAC, 47 – 63 Hz
PC interface	Ethernet, USB, RS 232

## Options

- Repetition rate 40 MHz (M40) or customized (some specifications may change)
  - Further two oscillator outputs for seeding purposes (on total 3 FC/APC outputs)
  - System without oscillator (AMP) for multiple beam extension systems (FC/APC input)
  - Variable laser repetition rate VAR (for more information see options page)
  - Phase-locked loop laser synchronization electronics LRC (see options page)
- Contact TOPTICA for customized systems


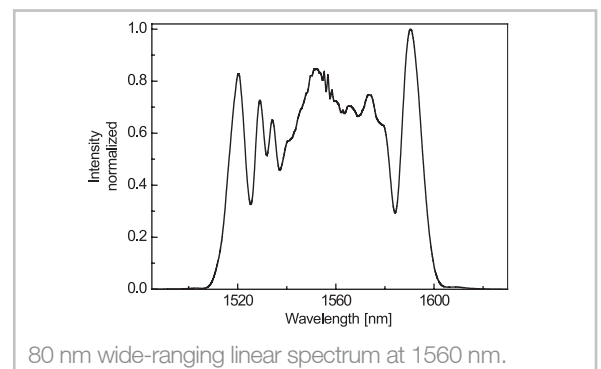
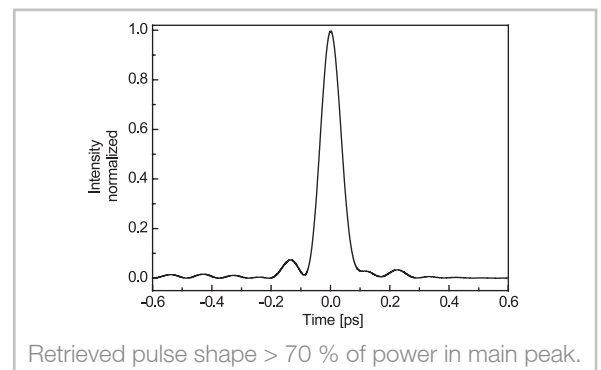
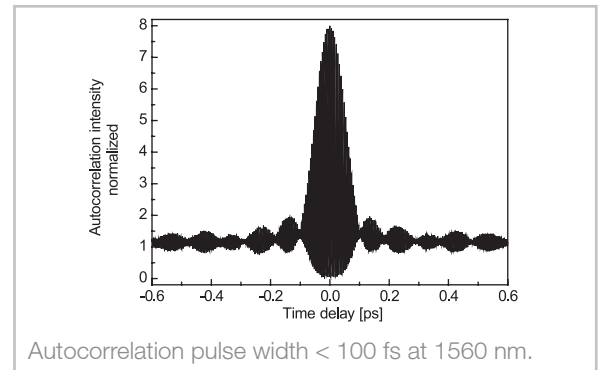


# FemtoFiber pro

**NIR**  
Near InfraRed

## Key features

- Fiber laser system with SAM mode-locked oscillator, core-pumped high power amplifier and second harmonic generation in a single box
- Control unit with laser drivers, power PC control, communication interfaces and power supply
- Motorized prism compressor for pulse optimization
- Free-beam output for 1560 nm and 780 nm, mechanically switchable
- Secondary oscillator fiber output (FC/APC connector)
- High bandwidth monitor output (SMA connector)

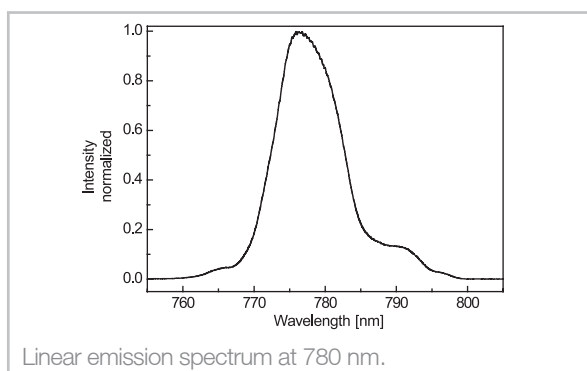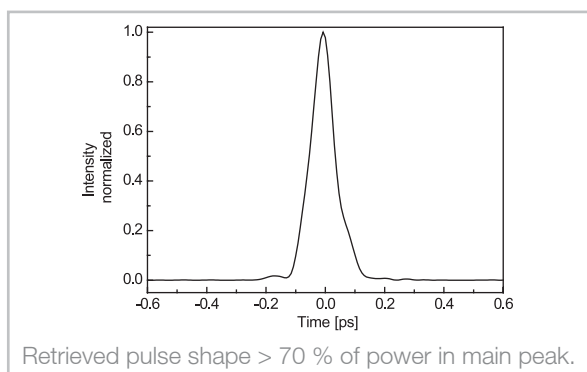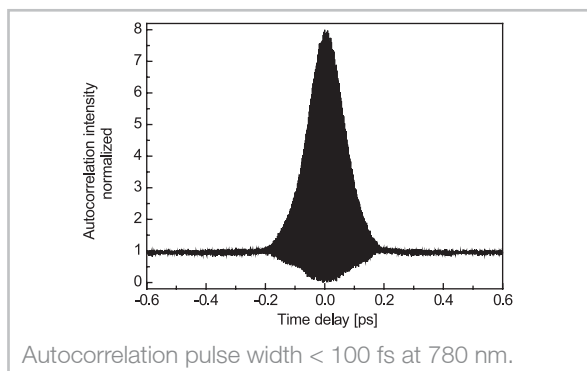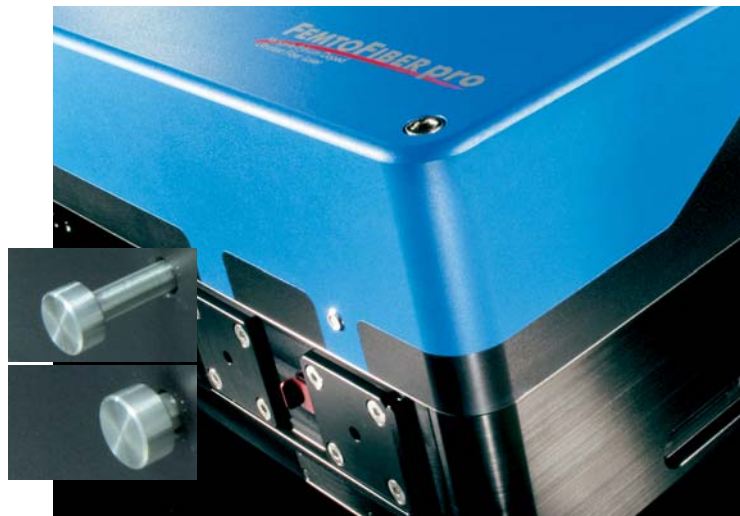

## Key specifications

Fundamental wavelength	<b>1560 nm</b>
Laser output power	> 350 mW
Pulse width	< 100 fs
Repetition rate	80 MHz standard
Second-harmonic wavelength	<b>780 nm</b>
Laser output power	> 140 mW
Linear polarization	> 95 % (horizontal)
Beam shape	TEM <sub>00</sub> , M <sup>2</sup> < 1.2
Beam size (1/e <sup>2</sup> )	Typ. Ø 1.2 mm (780 nm) Typ. Ø 3.5 mm (1560 nm)
Beam divergence	< 1 mrad (780 nm) < 2 mrad (1560 nm)
Output coupling	Free space
Laser head dimensions	280 x 229 x 151 mm <sup>3</sup> (w x d x h)
Laser head weight	< 10 kg
Control unit dimensions	235 x 315 x 140 mm <sup>3</sup> (w x d x h)
Control unit weight	< 4.5 kg
Line input	90 – 260 VAC, 47 – 63 Hz
PC interface	Ethernet, USB, RS 232

## Options

- Repetition rate 40 MHz (M40) or customized (some specifications may change)
  - Further two oscillator outputs for seeding purposes (on total 3 FC/APC outputs)
  - System without oscillator (AMP) for multiple beam extension systems (FC/APC input)
  - Variable laser repetition rate VAR (for more information see options page)
  - Phase-locked loop laser synchronization electronics LRC (see options page)
- Contact TOPTICA for customized systems



# FemtoFiber pro

## SCIR

### SuperContinuum InfraRed

#### Key features




- Fiber laser system with SAM mode-locked oscillator, core-pumped high power amplifier and supercontinuum generation in a single box
- Supercontinuum generated by highly nonlinear fiber
- Control unit with laser drivers, power PC control, communication interfaces and power supply
- Motorized prism compressor for pulse optimization
- Free-beam output
- Secondary oscillator fiber output (FC/APC connector)
- High bandwidth monitor output (SMA connector)
- Motorized chirp adjustment

#### Key specifications

Wavelength range continuum	<b>980 – 2200 nm</b>
Laser output power	> 150 mW
Pulse width	< 100 fs
Repetition rate	80 MHz standard
Linear polarization	> 95 % (horizontal)
Beam shape	TEM <sub>00</sub> , M <sup>2</sup> < 1.2
Beam size (1/e <sup>2</sup> )	Typ. Ø 4 mm
Beam divergence	< 1 mrad
Output coupling	Free space
Laser head dimensions	280 x 229 x 151 mm <sup>3</sup> (w x d x h)
Laser head weight	< 10 kg
Control unit dimensions	235 x 315 x 140 mm <sup>3</sup> (w x d x h)
Control unit weight	< 4.5 kg
Line input	90 – 260 VAC, 47 – 63 Hz
PC interface	Ethernet, USB, RS 232

#### Options

- Repetition rate 40 MHz (M40) or customized (some specifications may change)
  - Further two oscillator outputs for seeding purposes (on total 3 FC/APC outputs)
  - System without oscillator (AMP) for multiple beam extension systems (FC/APC input)
  - Variable laser repetition rate VAR (for more information see options page)
  - Phase-locked loop laser synchronization electronics LRC (see options page)
- Contact TOPTICA for customized systems

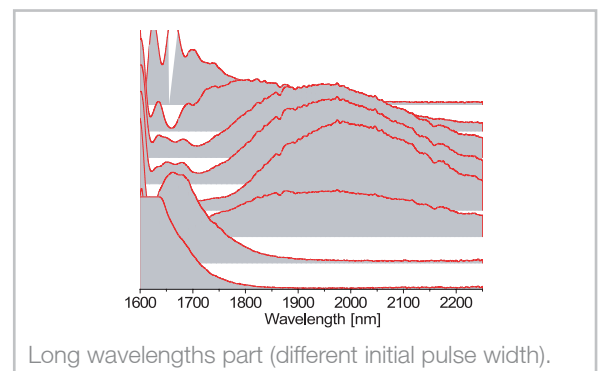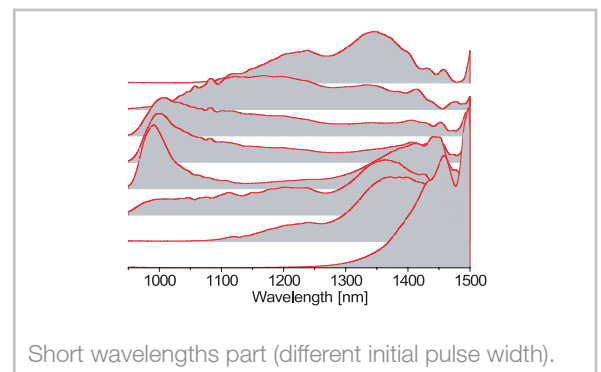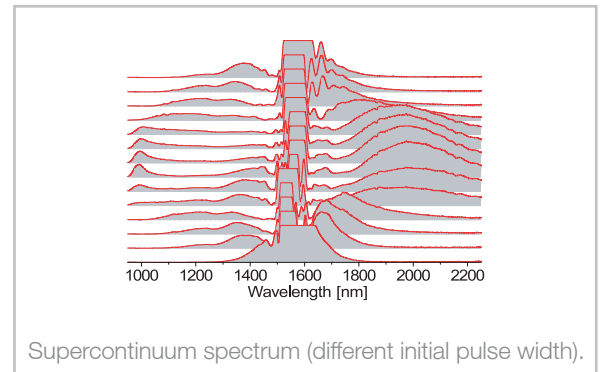





# FemtoFiber pro

## UCP

### Ultra Compressed Pulse

#### Key features

- Fiber laser system with SAM mode-locked oscillator and core-pumped high power amplifier
- Supercontinuum generated by highly nonlinear fiber
- Two motorized prism compressors for supercontinuum and pulse optimization
- Free-beam output
- Control unit with laser drivers, power PC control, communication interfaces and power supply
- Secondary oscillator fiber output (FC/APC connector)
- Interlock capabilities
- High bandwidth monitor output (SMA connector)

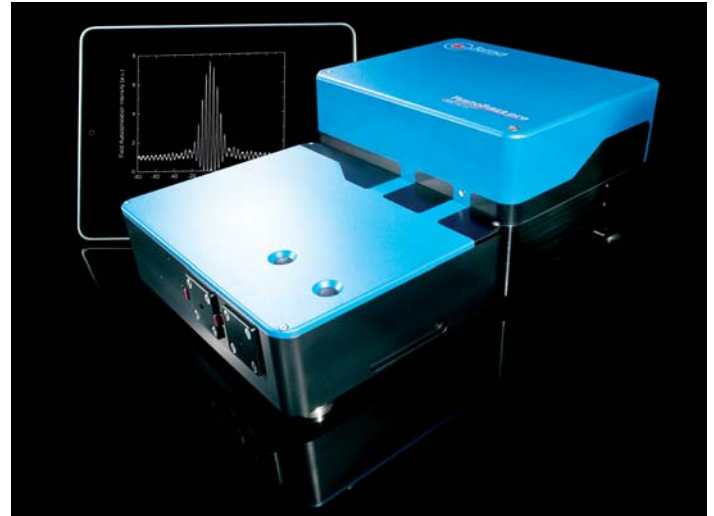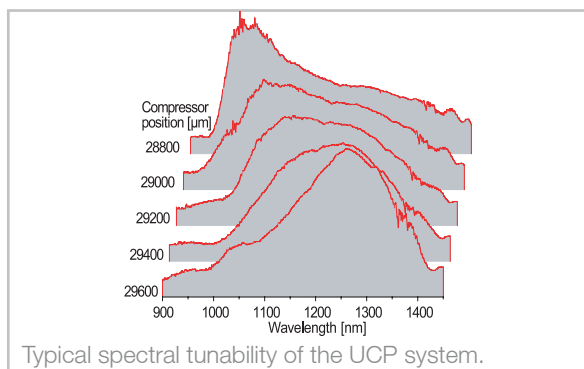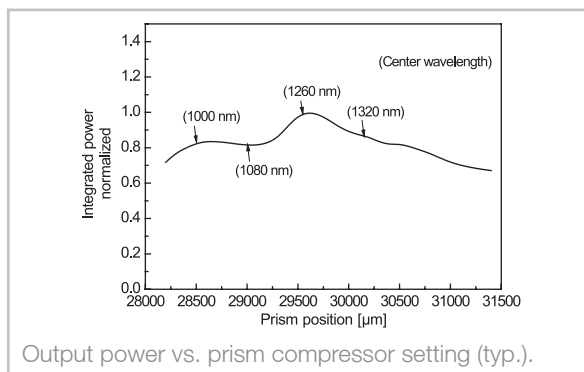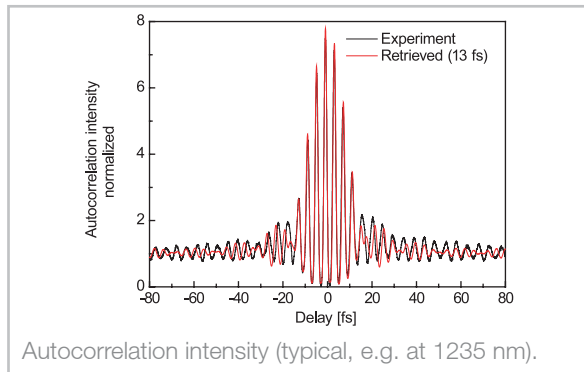

#### Key specifications

Wavelength range continuum	<b>980 – 1400 nm</b>
Laser output power	30 mW*
Pulse width	< 25 fs*
Repetition rate	80 MHz standard
Linear polarization	> 95 % (horizontal)
Beam shape	TEM <sub>00</sub> , M <sup>2</sup> < 1.2
Beam size (1/e <sup>2</sup> )	Typ. Ø 4 mm
Beam divergence	< 1 mrad
Output coupling	Free space
Laser head dimensions	280 x 458 x 151 mm <sup>3</sup> (w x d x h)
Laser head weight	< 15 kg
Control unit dimensions	235 x 315 x 140 mm <sup>3</sup> (w x d x h)
Control unit weight	< 4.5 kg
Line input	90 – 260 VAC, 47 – 63 Hz
PC interface	Ethernet, USB, RS 232
*) Dependent on center wavelength	

#### Options

- Repetition rate 40 MHz (M40) or customized (some specifications may change)
- Further two oscillator outputs for seeding purposes (on total 3 FC/APC outputs)
- System without oscillator (AMP) for multiple beam extension systems (FC/APC input)
- Can be combined with TNIR option (contact TOPTICA)
- Variable laser repetition rate VAR (for more information see options page)
- Phase-locked loop laser synchronization electronics LRC (see options page)

Contact TOPTICA for customized systems



# FemtoFiber pro

## TVIS

Tunable VISible

### Key features

- Fiber laser system with SAM mode-locked oscillator and core-pumped high power amplifier
- Supercontinuum generated by highly nonlinear fiber
- Two motorized prism compressors for supercontinuum and pulse optimization
- Second harmonic generation by manually tunable crystal
- Free-beam output
- Control unit with laser drivers, power PC control, communication interfaces and power supply
- Secondary oscillator fiber output (FC/APC connector)
- High bandwidth monitor output (SMA connector)




### Key specifications

Wavelength range continuum	<b>488 – 640 nm</b>
Laser output power	1 – 10 mW
Pulse width	< 1 ps
Repetition rate	80 MHz standard
Linear polarization	> 95 % (horizontal)
Beam shape	TEM <sub>00</sub> , M <sup>2</sup> < 1.2
Beam size (1/e <sup>2</sup> )	t.b.d.
Beam divergence	< 1 mrad
Output coupling	Free space
Laser head dimensions	280 x 458 x 151 mm <sup>3</sup> (w x d x h)
Laser head weight	< 15 kg
Control unit dimensions	235 x 315 x 140 mm <sup>3</sup> (w x d x h)
Control unit weight	< 4.5 kg
Line input	90 – 260 VAC, 47 – 63 Hz
PC interface	Ethernet, USB, RS 232

### Options

- Repetition rate 40 MHz (M40) or customized (some specifications may change)
- Further two oscillator outputs for seeding purposes (on total 3 FC/APC outputs)
- System without oscillator (AMP) for multiple beam extension systems (FC/APC input)
- Can be combined with TNIR option (contact TOPTICA)
- Variable laser repetition rate VAR (for more information see options page)
- Phase-locked loop laser synchronization electronics LRC (see options page)

Contact TOPTICA for customized systems

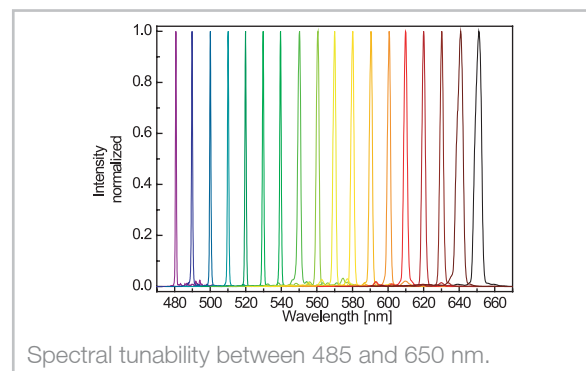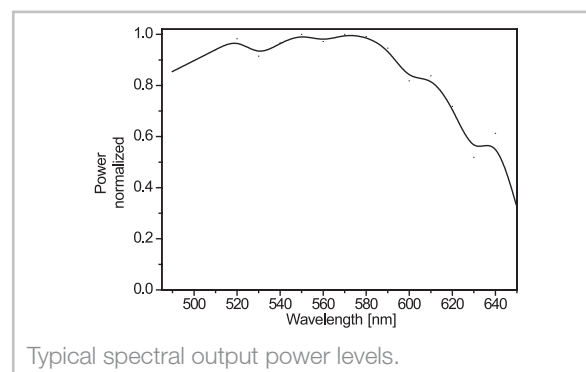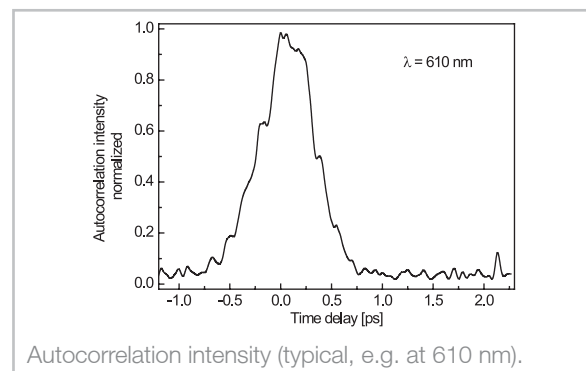







# FemtoFiber pro

## TNIR

Tunable Near InfraRed

### Key features

- Fiber laser system with SAM mode-locked oscillator and core-pumped high power amplifier
- Supercontinuum generated by highly nonlinear fiber
- Motorized prism compressor for pulse optimization
- Second harmonic generation by manually tunable crystal
- Free-beam output
- Control unit with laser drivers, power PC control, communication interfaces and power supply
- Secondary oscillator fiber output (FC/APC connector)
- High bandwidth monitor output (SMA connector)

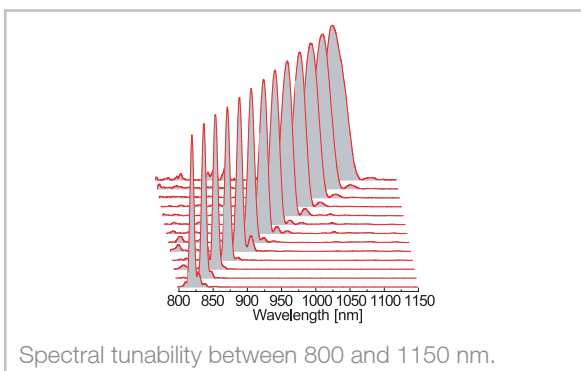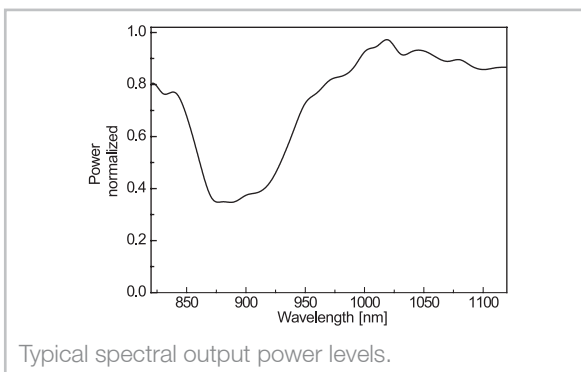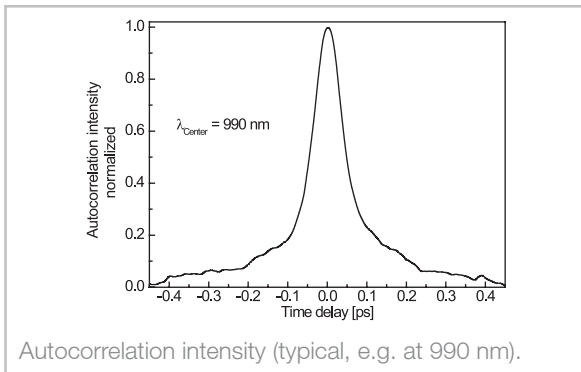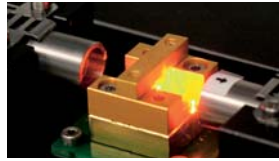

### Key specifications

Wavelength range continuum	<b>830 – 1100 nm</b>
Laser output power	1 – 5 mW
Pulse width	< 200 fs
Repetition rate	80 MHz standard
Linear polarization	> 95 % (horizontal)
Beam shape	TEM <sub>00</sub> , M <sup>2</sup> < 1.2
Beam size (1/e <sup>2</sup> )	t.b.d.
Beam divergence	< 1 mrad
Output coupling	Free space
Laser head dimensions	280 x 458 x 151 mm <sup>3</sup> (w x d x h)
Laser head weight	< 15 kg
Control unit dimensions	235 x 315 x 140 mm <sup>3</sup> (w x d x h)
Control unit weight	< 4.5 kg
Line input	90 – 260 VAC, 47 – 63 Hz
PC interface	Ethernet, USB, RS 232

### Options

- Repetition rate 40 MHz (M40) or customized (some specifications may change)
- Further two oscillator outputs for seeding purposes (on total 3 FC/APC outputs)
- System without oscillator (AMP) for multiple beam extension systems (FC/APC input)
- Can be combined with TVIS or UCP option (contact TOPTICA)
- Variable laser repetition rate VAR (for more information see options page)
- Phase-locked loop laser synchronization electronics LRC (see options page)

Contact TOPTICA for customized systems

# FemtoFiber pro

## Options

### Option: VAR — VArIable laser Repetition rate

- Adaptation to the oscillator unit, enabling modulation of the repetition rate
- Adjustable resonator length by Piezo transducer and translation stage
- Piezo resonance frequency > 5 kHz
- Repetition rate tuning range > 200 kHz

### Option: LRC — LAser Repetition rate Control

- Phase-locked loop electronics for synchronization of the laser pulse train to external reference signal
- RMS jitter t.b.d.
- Electronics rack with power supply
- USB interface and control software
- ECOPS capability

### Option: M40 — 40 MHz repetition rate


- Oscillator design with 40 MHz repetition rate

Some specifications may change


### Option: AMP — MUlti beam configuration

- System without oscillator for multi beam systems
- FC/APC fiber input


## Electronically controlled optical sampling — ECOPS

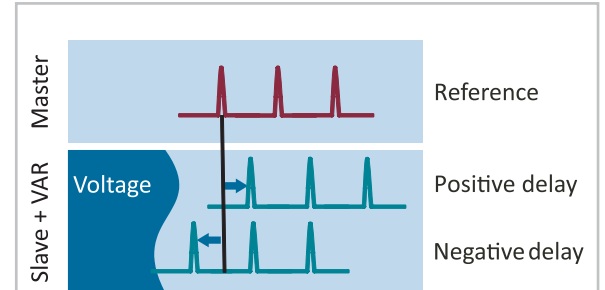

ECOPS: Mode of operation.

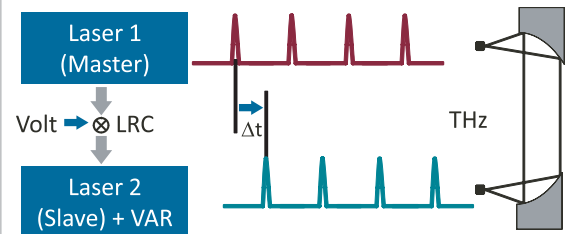


ECOPS: Example implementation for THz.

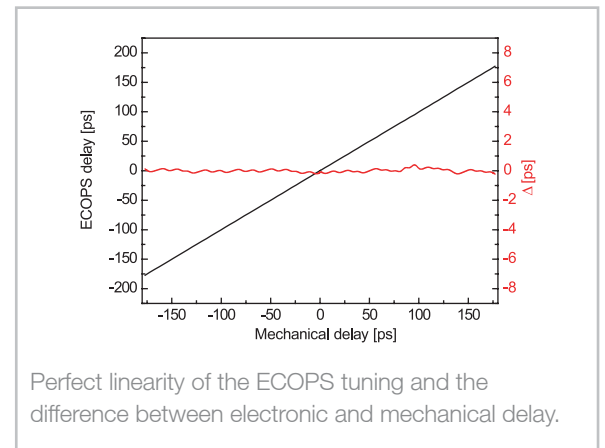

Perfect linearity of the ECOPS tuning and the difference between electronic and mechanical delay.

## Fiber coupling to other units

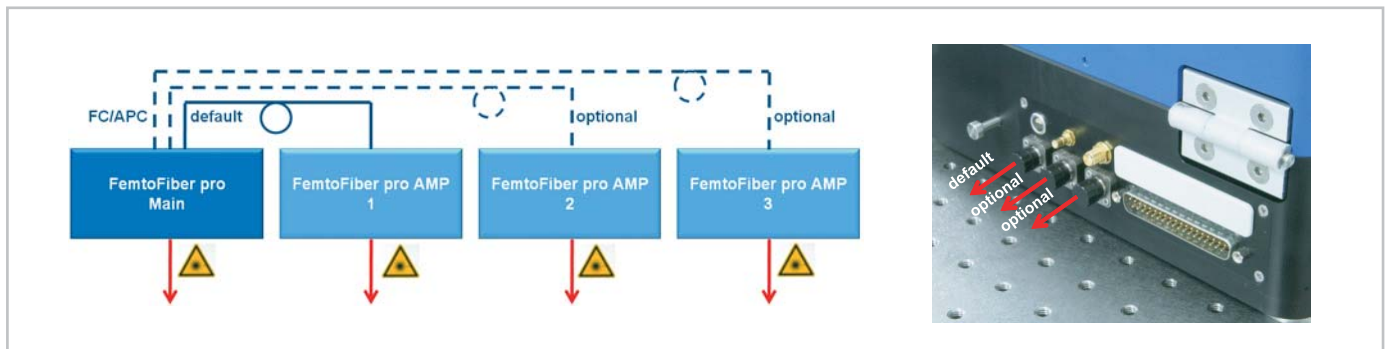

# FemtoFiber pro Control Unit

## Key features


- Built-in power PC for system control
- Easy communication through web browser
- Access to motorized controls, such as variable pulse compression
- LabVIEW™ routines available for system integration
- Manual interface: Push ON/OFF button only
- Key lock switch
- Interlock capabilities
- 12 inch rack housing including interfaces, driver electronics for pump diodes and power supplies




FemtoFiber pro control unit.

The laser control unit includes a power PC that controls all laser parameters. This ensures turnkey and hands-off operation with a single ON/OFF button for the user. Standard communication interfaces (Ethernet, USB and RS 232) give access to all relevant parameters for an easy integration in complex setups.


The user can choose between three alternatives. Simplest way is to login via any web browser already installed on the user's computer. Another way is employing the LabVIEW™ routine that is included free of charge. The most sophisticated way is to use self-written scripts with a pool of pre-defined commands.

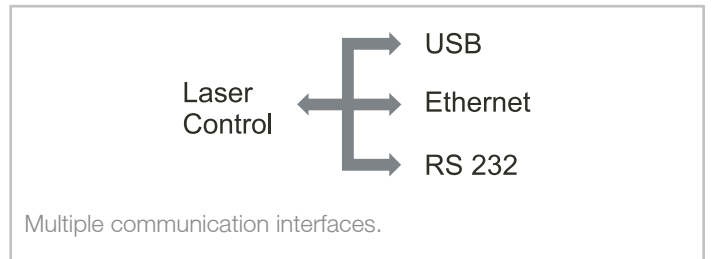

Multiple communication interfaces.

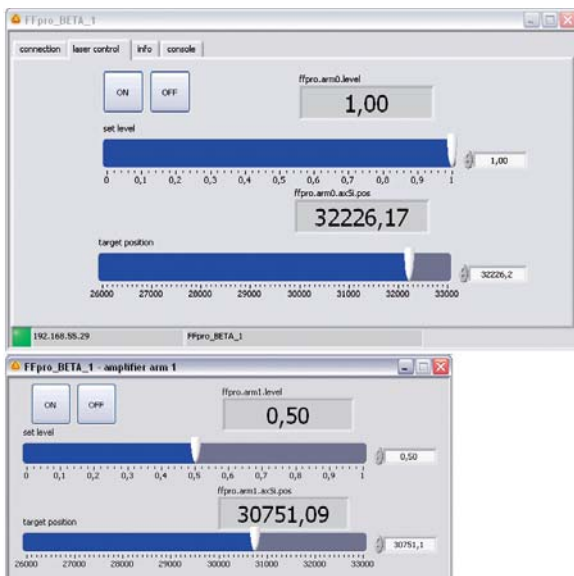

Graphical user interface based on LabVIEW™.

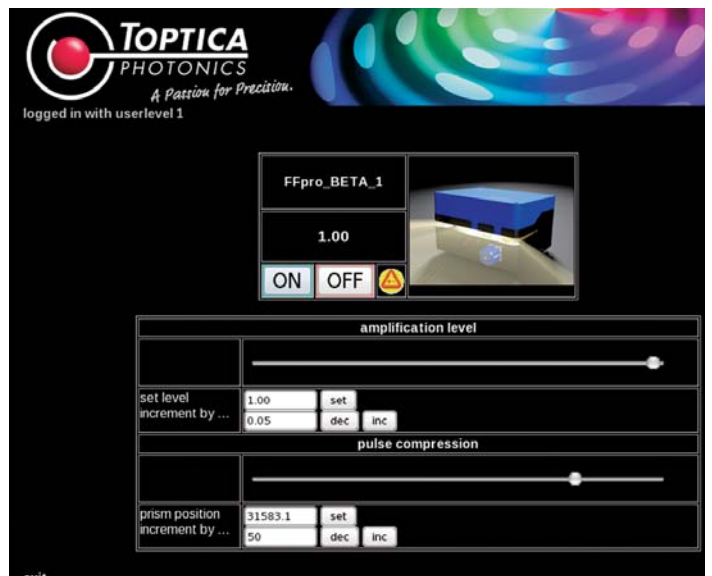


Web browser based GUI.

# All-Color Fiber Laser iChrome



## Key specifications

Wavelength range	<b>488 – 640 nm</b>
Average fiber output power	> 1.5 mW
Repetition rate	40 MHz ± 1 MHz
Bandwidth (FWHM)	< 3 nm
Pulse width	3.5 ps ± 1 ps
Polarization	Linear, 100 : 1
Trigger output	Typ. - 300 mV, 1 ns
Spatial mode	TEM <sub>00</sub>
Fiber delivery	PM SM fiber, 2 m length MFD 4 μm @ 480 nm
Accuracy	± 1 nm
Speed	≥ 50 nm/s
Line-up time	≤ 500 ms
Warm-up time	2 h
Operating temperature	22°C ± 3°C
Air humidity	< 70 % non-condensing
Power consumption	< 75 W
Dimensions	450 x 570 x 230 mm <sup>3</sup> (w x d x h)
Weight	21.5 kg

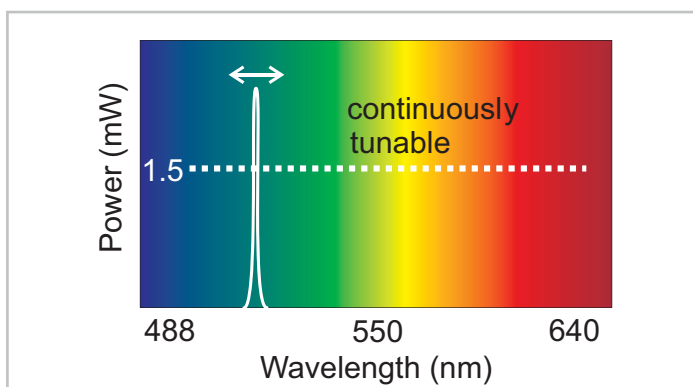

## Applications

- Pump probe spectroscopy
- Fluorescence life time measurement
- Excitation spectra measurement
- Confocal microscopy

## Flexible fiber laser

The newly introduced iChrome is a fiber laser with the flexibility to set automatically the laser output to any wavelength in the visible – from 488 nm to 640 nm. In contrast to conventional white light sources, the narrow bandwidth laser pulses are not filtered out from an intrinsically noisy supercontinuum. Moreover, even the coherence of the fundamental laser is preserved during the frequency conversion processes. This ensures that the visible light exhibits the best intensity noise performance. Additionally all optical components are polarization-maintaining. This results in a stable linear polarization of the output beam.

## Fiber coupled output

The laser output of the iChrome is delivered by a single-mode and polarization-maintaining fiber. Independent from the chosen wavelength the fiber output exhibits a smooth TEM<sub>00</sub> profile, an excellent beam product with  $M^2 < 1.1$  and a linear polarization with a PER\* > 1:100. The fiber makes the connection between the high quality light source and the experiment as simple as possible. Pointing stability is guaranteed and all colors are delivered by the same fiber.

## Fully automated operation

The entire laser system is extremely user friendly: No alignment procedures of any optical components distract the user from the main task – to produce results. Besides the fiber delivery also the full automation of the premium light source helps to avoid wasting time. The built-in power PC controls all necessary parameters to ensure smooth operation every minute, hour, day, week and month. The power PC also hosts a webserver and is equipped with an ethernet connection. Therefore, all user commands can be sent from an ordinary web browser to the laser system.

\*PER = Polarization Extinction Ratio

# Industrial-Grade Fiber Lasers

## PicoFYb 1030 nm, 1064 nm / FemtoFErb 1560 nm

### Applications

- Micro material processing
- Time-domain Terahertz generation
- Light source for microscopy
- Medical surgery / examination
- Ophthalmology

### Ultrafast, ultra stable, turnkey

The PicoFYb / FemtoFErb laser systems are turnkey fiber-based picosecond and femtosecond mode-locked fiber oscillators. The PicoFYb is ideally suited for seeding industrial laser systems. The FemtoFErb is the most compact and cost-effective source for Terahertz generation.



### Key features

The PicoFYb / FemtoFErb laser pulses exhibit excellent amplitude and frequency jitter parameters. The PicoFYb is typically amplified to multi-Watt levels in the MOPA (master oscillator, power amplifier) laser systems of our customers.

The FemtoFErb has an output power of more than 100 mW and fits perfect to applications such as Terahertz generation or metrology.

Based on state-of-the-art FemtoFiber technology, the PicoFYb / FemtoFErb lasers set the benchmark for the new

generation of ultrafast lasers satisfying the most advanced industrial requirements.


Applications benefit from the most stable, compact and cost effective laser design, but in particular from the unique fiber output coupling, a key to a modular customer integration.

The specifications are derived from our current OEM developments. Further customization (wavelength, pulse duration, footprint, interfacing, etc.) is possible under usual OEM conditions.


### Key specifications

Model	PicoFYb		FemtoFErb
	1030 nm AMP	1064 nm AMP	1560 nm AMP
Wavelength	1030.5 ± 0.5 nm	1064.3 ± 0.5 nm	1560
Average power	> 10 mW <sup>*/**</sup>	10 mW <sup>*/**</sup>	> 100 mW <sup>*/**</sup>
Repetition rate	20 MHz ± 0.5 MHz	20 MHz ± 0.5 MHz	50 MHz / 100 MHz
Pulse duration (FWHM)	< 10 ps <sup>**</sup>	6 ± 2 ps <sup>**</sup>	< 150 fs <sup>**</sup>
Spectral width	< 0.5 (0.3) nm	< 0.5 (0.3) nm	< 40 nm
TBP	Typ. 0.5	Typ. 0.5	tbd
RF side mode suppression	< -60 dBC	< -60 dBC	tbd
P-p-noise	< 3 %	< 3 %	tbd
Output	SM PM fiber (narrow key FC/APC receptacle)		
Dimensions	122 x 202 x 69 mm <sup>3</sup> (w x d x h)		


<sup>\*</sup>) Fixed output power, not adjustable    <sup>\*\*</sup>) After fiber delivery

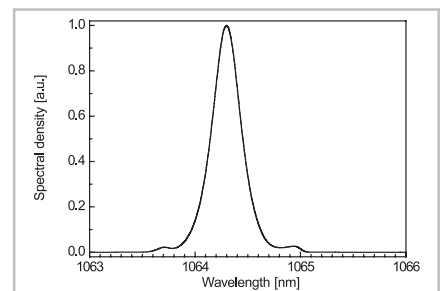

PicoFYb 1064 AMP: Linear spectrum with less than 0.5 nm linewidth.

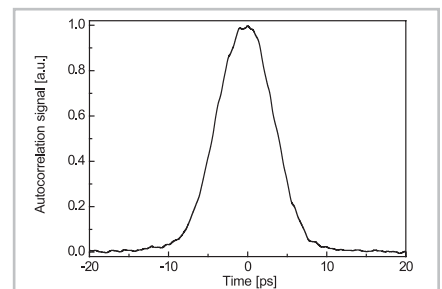


PicoFYb 1064 AMP: Autocorrelation data with less than 10 ps pulse width.

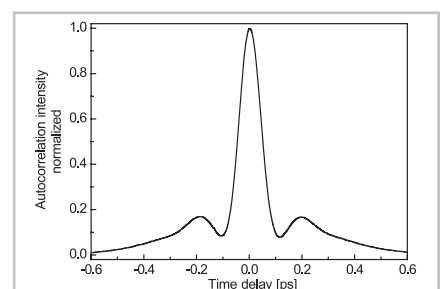

FemtoFErb 1560 AMP: Autocorrelation data.

# Distributors

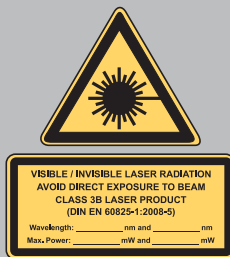

Australia & New Zealand  
**Lastek Pty. Ltd.**  
Mr. Alex Stanco  
Thebarton Campus  
University of Adelaide  
10 Reid Street  
5031 Thebarton, SA  
Australia  
Phone: +61 8 8443 8668  
Fax: +61 8 8443 8427  
sales@lastek.com.au  
www.lastek.com.au

China  
**Universal (Hong Kong)  
Technology Co. Ltd.**  
Mr. Alex Cai  
Room 615-619  
Capital Group Plaza  
No. 6 Chaoyangmen Beidajie  
Beijing 100027, P.R. China  
Phone: +86 10 8528 3377  
Fax: +86 10 8528 3344  
oe@universalkco.com.cn  
www.universalkco.com.cn

France  
**Opton Laser International**  
Dr. Costel Subran  
Parc Club d'Orsay Université  
29, rue Jean Rostand  
F-91893 Orsay Cedex, France  
Phone: +33 1 6941 0405  
Fax: +33 1 6941 3290  
ventes@optonlaser.com  
www.optonlaser.com

India  
**Simco Global Technology &  
Systems Ltd.**  
Dr. R. S. Daryan  
Simco House (Head Office)  
14 Bhawani Kunj  
Behind Sector D-II Vasant Kunj  
110017 New Delhi, India  
Phone: +91 11 2689 9867  
Fax: +91 11 2612 4461  
simcorsd@del2.vsnl.net.in  
www.simco-groups.com

Israel  
**Lahat Technologies Ltd.**  
Mr. Kfir Ben-Yehuda  
Teradion Industrial Zone  
M.P. Misgav, 20179, Israel  
Phone: +972 4 999 0151  
Fax: +972 4 999 0826  
sales@lahat.co.il  
www.lahat.co.il

Japan  
**INDECO, INC.**  
Mr. Kazunobu Takahashi  
1-11-14, Kasuga, Bunkyo-ku  
Tokyo 112-0003, Japan  
Phone: +81 3 3818 4011  
Fax: +81 3 3818 4015

**Sonezaki Sawada Bld. 6F**  
2-1-13, Sonezakishinchi,  
Kita-ku,  
Osaka 530-0002, Japan  
Phone: +81 6 6341 5799  
Fax: +81 6 6341 5798

ind@indecop.jp  
www.indecop.jp

Korea  
**JINSUNG LASER**  
Mr. Ha-Won Lee  
#535-5, Bongmyung-Dong  
Yusung-Gu  
Hanjin Officetel Rm# 1016  
Daejeon, 305-301, South Korea  
Phone: +82 42 823 5300  
Fax: +82 42 823 7447  
sales@jinsunglaser.com  
www.jinsunglaser.com

Singapore / Malaysia  
**Precision Technologies Pte Ltd**  
Mr. Vincent Chan  
211 Henderson Road # 13-02  
Henderson Industrial Park  
Singapore 159552  
Phone: +65 6273 4573  
Fax: +65 6273 8898  
precision@pretech.com.sg  
www.pretech.com.sg

Taiwan  
**SLEO Photonics Co. Ltd.**  
Mr. Jimmy Chao  
6F, No. 2, Lane 74  
An-der Street  
Hsin Tien City, Taipei County  
Taiwan 231, R.O.C.  
Phone: +886 2 2211 5408  
Fax: +886 2 2211 5401  
sleo.jimmy@msa.hinet.net

Turkey  
**CS Analytical Laboratory  
Equipments Limited**  
Mr. Ercan Acikalın  
Ergin Sokak No: 27/5  
06580 Mebusevler - Tandoğan  
Ankara, Turkey  
Phone: +90 312 223 0302  
Fax: +90 312 223 0305  
can@csanalitik.com.tr  
www.csanalitik.com.tr

United Kingdom & Ireland  
**Mr. Howard Potter**  
Unit 4H Lansbury Estate  
Woking, Surrey, UK, GU21 2EP  
Great Britain  
Phone: +44 1483 799 030  
Fax: +44 1483 799 076  
howard.potter@toptica.com  
www.toptica.com

USA & Canada  
**TOPTICA Photonics, Inc.**  
1286 Blossom Drive  
Victor/Rochester, NY 14564  
U.S.A.  
Phone: +1 585 657 6663  
Fax: +1 877 277 9897  
sales@toptica.com  
www.toptica.com

Every other country  
not listed:

**TOPTICA Photonics AG**  
Lochhamer Schlag 19  
D-82166 Graefelfing / Munich  
Germany  
Phone: +49 89 85837-0  
Fax: +49 89 85837-200  
sales@toptica.com  
www.toptica.com



**TOPTICA Photonics AG**  
Lochhamer Schlag 19  
D-82166 Graefelfing / Munich  
Germany  
Phone: +49 89 85837-0  
Fax: +49 89 85837-200  
info@toptica.com  
www.toptica.com

**TOPTICA Photonics, Inc.**  
1286 Blossom Drive  
Victor/Rochester, NY 14564  
U.S.A.  
Phone: +1 585 657 6663  
Fax: +1 877 277 9897  
sales@toptica.com  
www.toptica.com