



## ARI Gas Chromatograph

*A field-deployable modular GC for time-of-flight mass spectrometry (TOF-MS)*


### Features:

**Flexibility** – up to three GC columns depending on application, pre-concentration methods (e.g. thermal desorption vs cryo-mechanical), detector type(s)

**Field deployable** – ruggedized GC with small, light-weight, footprint for automated operation

**Detector compatibility** – true plug-and-play operation with all ARI TOF-MS systems (hardware and software, including Tofwerk Acquity)

**Dual detectors** – designed to operate automatically with two detectors for a more complete data set

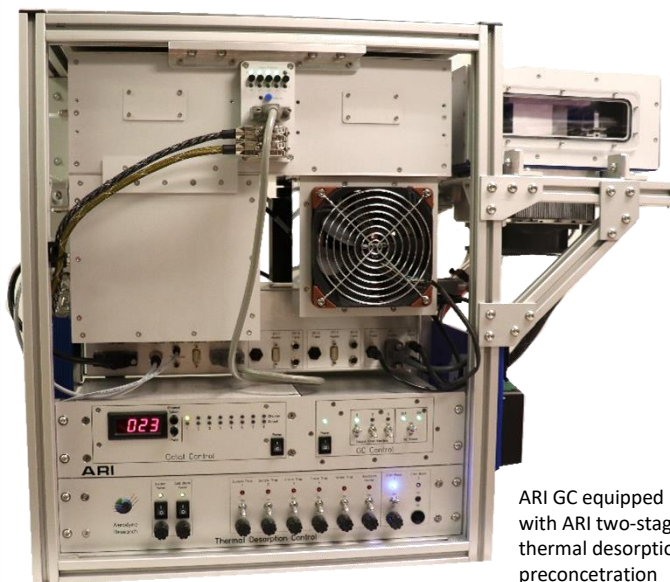

ARI GC equipped with ARI two-stage thermal desorption preconcentration system

**Turn-Key** – ready-to-run on delivery due to pre-installed analytical column(s) with custom pre-loaded method(s)

**Speed** – flow path allows user to perform multiple GC operations at once; the benefit of this parallel operation is a faster GC cycle (10 – 30 min) without loss of performance

### GC Technology

Aerodyne Research, Inc. (ARI) offers a modular gas chromatograph (GC) that combines recent innovations in fast GC separation with highly selective and sensitive detection in a field-deployable package for the measurement of VOCs and OVOCs.

Aerodyne works with our customers to provide a GC system with a **custom separation method** appropriate for their analytical needs. The instrument delivers with the required capillary column installed and analytical methods pre-loaded into the Windows-based control software, for a true **turn-key** operation.

### Data Analysis Software

Wavemetrics Igor Pro analysis software (TERN) allows for retention time correction, automated baseline and peak fitting. TERN provides de-convoluted peak fitting for overlapping chromatographic peaks. Compatible for analysis of both high resolution (HR) and unit mass resolution (UMR) data.

### Compatible Accessories

#### Sample Preconcentrators

- ARI Thermal Desorption (TDPC)
- ARI Cryomechanical (CMPC)
- Other commercial preconcentration systems


#### Direct Injection Methods

- Passivated sample loop (0.005 – 1 mL)
- Split/Splitless Injection Port
- Commercial autosamplers

#### Detectors

- Vocus PTR-TOF-MS
- EI-TOF-MS
- CIMS-TOF-MS
- Other commercial detectors (Q-MS, FID, ECD, etc.)

## Performance Specifications:

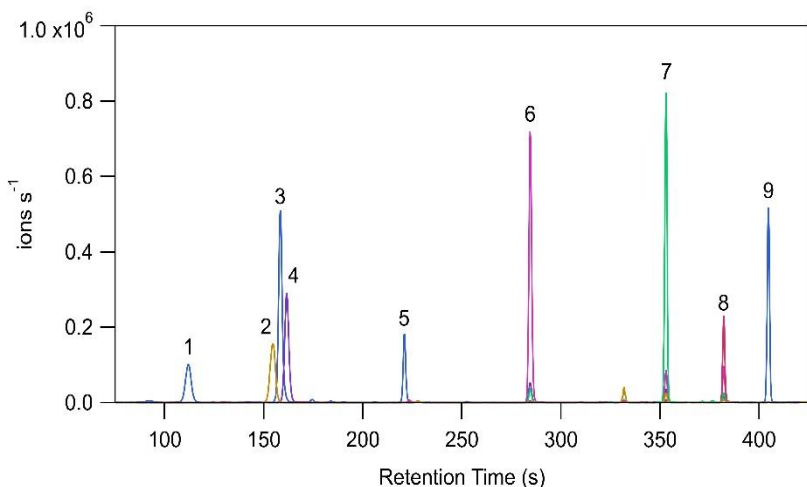

### Column Oven Specifications

Max heating rate: 120 °C/min

Max cooling rate: 250 °C/min

Max temperature: 260 °C

Temperature accuracy: ± 0.5 °C (ramped);  
± 0.25 °C (isothermal)

Column compatibility: all fused silica and metal columns  
(60m @ ≤0.32mm ID, 30m @ ≥0.53mm ID)

Capillary columns are swappable without venting  
detector

### Sampling / Reagent Gas Specifications

3-inlet system for analysis of sample (ambient), zero, and  
calibration gases

Sample gas: -40 to 50 °C; 0-90% RH; 50-125 kPa

Carrier gas compatibility: helium, nitrogen, hydrogen

Sample flow rate: up to 200 sccm, ± 1% accuracy

Carrier gas flow rate: up to 10 sccm

Complete inert/passivated sample flow path

(e.g. PFA, SilcoTek Siltek/Sulfinert passivated stainless steel)

### Data Outputs

USB

### Weight, Size, Power

Weight: 24 kg

Dimensions: 55 cm x 55 cm x 30 cm

Max power: 600 W, 120/240 V, 50/60 Hz (start-up)  
300 W (typical operation)

Figure 1. ARI GC equipped with ARI two-stage thermal desorption preconcentration system with EI-TOF detection (TD-GC-EI-TOF). Analysis of multi-component calibration mixture in nitrogen. 1000 mL sample of 0.75 ppb each (1) n-hexane, (2) carbon tetrachloride, (3) 2,2,4-trimethylpentane, (4) benzene, (5) n-octane, (6) ethyl benzene, (7) 1,3,5-trimethylbenzene, (8) limonene, (9) n-undecane.

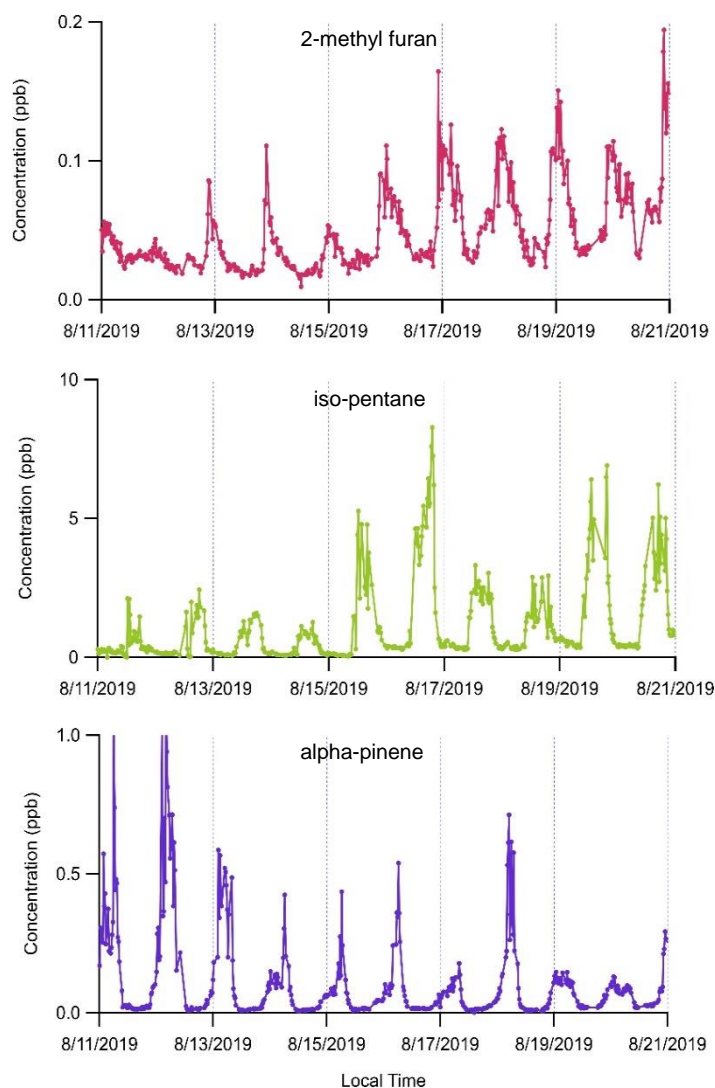


Figure 2. Ambient data acquired with ARI GC equipped with ARI two-stage thermal desorption preconcentration system with EI-TOF detection (TD-GC-EI-TOF)

*Aerodyne specializes in collaboration and custom design. Please contact us if you would like to discuss additional measurement options and applications.*

### REFERENCES:

Clafin, M. S., et al. (2021) Atmos. Meas. Tech. 14: 133-152. doi: 10.5194/amt-14-133-2021

Lerner, B. M., et al. (2017) Atmos. Meas. Tech. 10(1): 291-313. doi: 10.5194/amt-10-291-2017

Isaacman-Van Wertz et al. (2017) J Chromatogr A; 1529:81-92. doi: 10.1016/j.chroma.2017.11.005




## ARI TDPC-GC for TOF-MS

*A field-deployable, in situ thermal desorption preconcentration (TDPC) system for the Aerodyne gas chromatograph (GC) with time-of-flight mass spectrometry (TOF-MS).*

Mating the Aerodyne two-stage TDPC with the ARI GC expands the use of the GC system to a wide range of environmental applications. The system allows the user to automatically collect and analyze ambient air samples with sub-ppt sensitivity.


The first stage allows fast collection (maximum 150 sccm) of large samples (1 L of air, typical) with post-collection water purge of the sample in humid environments. The sample is then transferred to the second stage to focus the sample into a narrow injection onto the GC column.

Using a three 2-port chromatography valve flow path gives the user 8-modes of operation. These 8 modes allow various functions of the GC to be conducted in parallel without interfering in the analysis. For example, the user is able to collect sample onto the first stage of the TDPC while simultaneously analyzing a chromatogram.

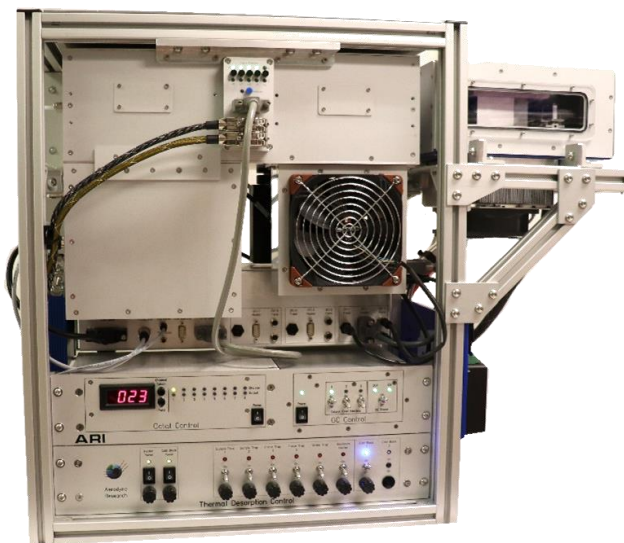

ARI GC equipped with ARI two-stage thermal desorption preconcentration system

*The benefit of this parallel operation is a faster GC cycle, 15 – 30 minutes, without the loss of performance.*

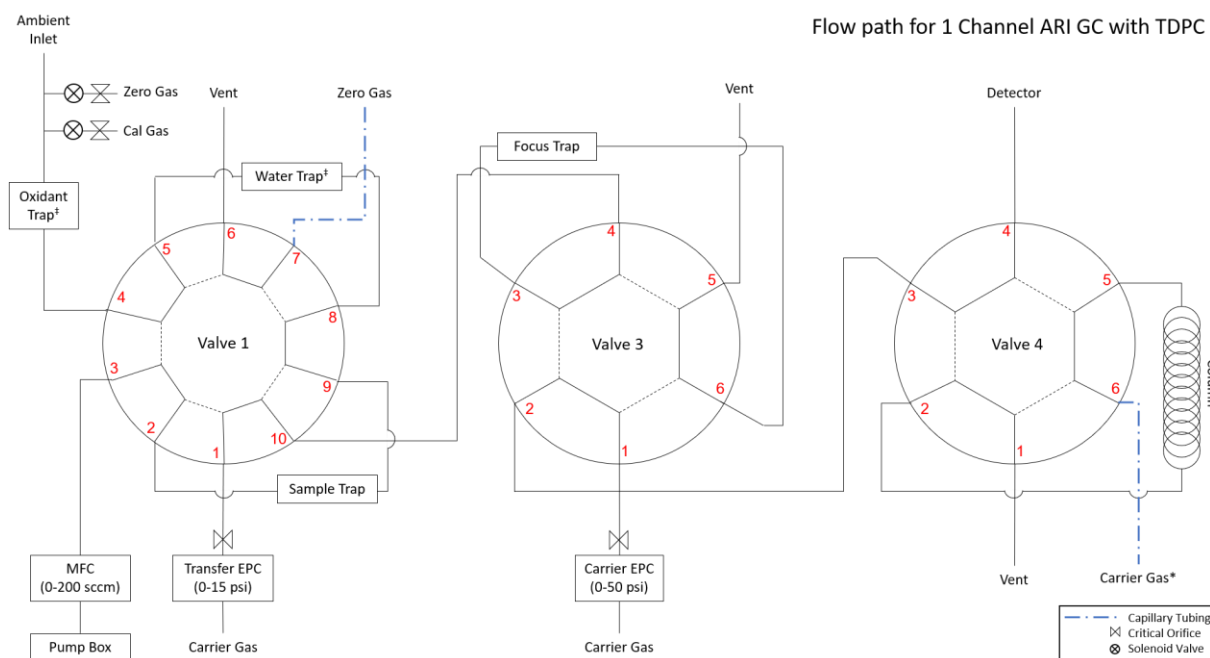


### *Sub parts-per-trillion sensitivity*

The coupling of the TDPC for online sampling with the chromatographic separation of the ARI GC results in limits of detection < 1 ppt for a wide range of atmospherically relevant compounds (see Appendix I) with average R<sup>2</sup> of 0.996 (see Appendix II).



## ARI TDPC-GC with 2 Column Separation

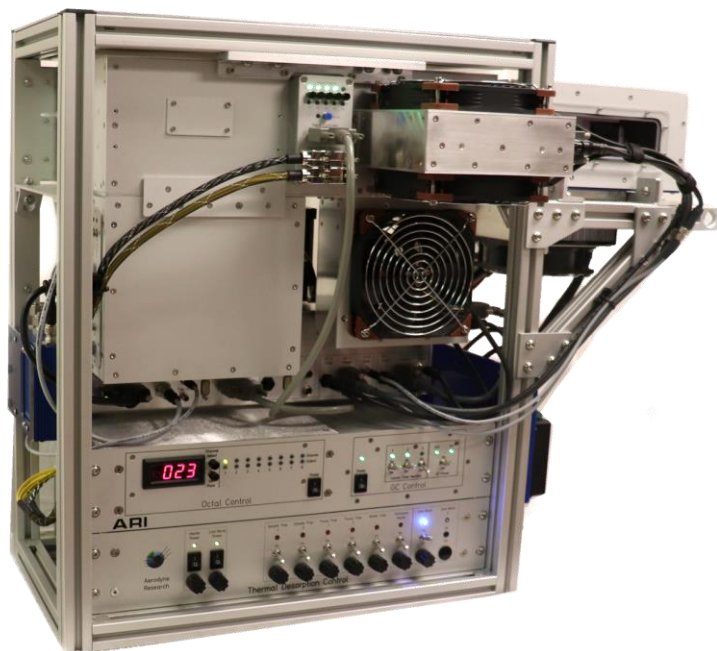

2-Channel ARI GC equipped with ARI two-stage thermal desorption preconcentration system

The ARI TDPC-GC for TOF-MS is available as a 2-channel system for the in situ analysis of VOCs and OVOCS.


The 2-channel system operates with a 2 column separation for the analysis of a wider volatility range.

Each GC separation channel is optimized for a different volatility range to expand the types of compounds resolved by the GC, while maintaining a fast GC cycle.

The 2-Channel ARI TDPC-GC system operates on a 30 minute cycle (typical), which includes sample preconcentration, post collection water purge, and a 10 minute chromatogram from each separation channel.

The user can choose a combination of volatility ranges that fits their analysis needs. The combination of Channels A and B will resolve  $C_2 - C_{12}$  VOCs, while the combination of Channels B and C will resolve  $C_5 - C_{20}$  VOCs.

If the users analysis needs change, either separation channel is easily modified by swapping in a different column assembly. This switch takes less than 15 minutes, and is done without venting the detector.

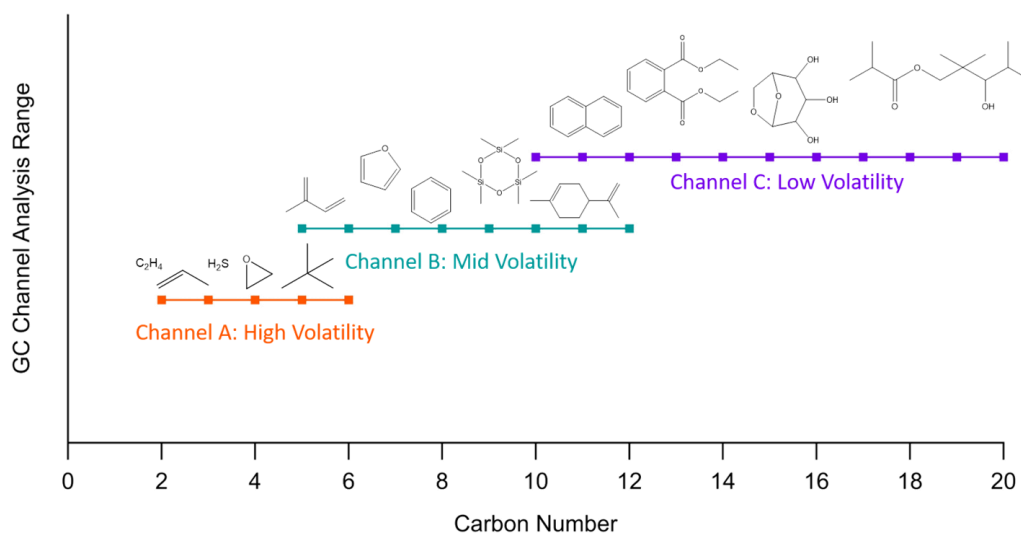

Example volatility ranges offered with the ARI GC. The 1-Channel ARI GC system is able to resolve one of these separation ranges, while the 2-Channel system can resolve a combination of Channels A, B, or C for the analysis of a wider range of VOCs and OVOCS.


The molecular structures shown are examples of compounds in the analysis range of each GC channel, but compounds that can be resolved by each channel are not limited to those shown.



## Exemplary ARI TDPC-GC-Vocus Chromatogram

Aerodyne TDPC-GC chromatographic separation of a selection of atmospherically relevant VOCs, detected with Vocus PTR-TOF-MS (proton transfer reaction time-of-flight mass spectrometry).

Figure shows a 10 minute chromatographic separation typically performed by the ARI GC system.

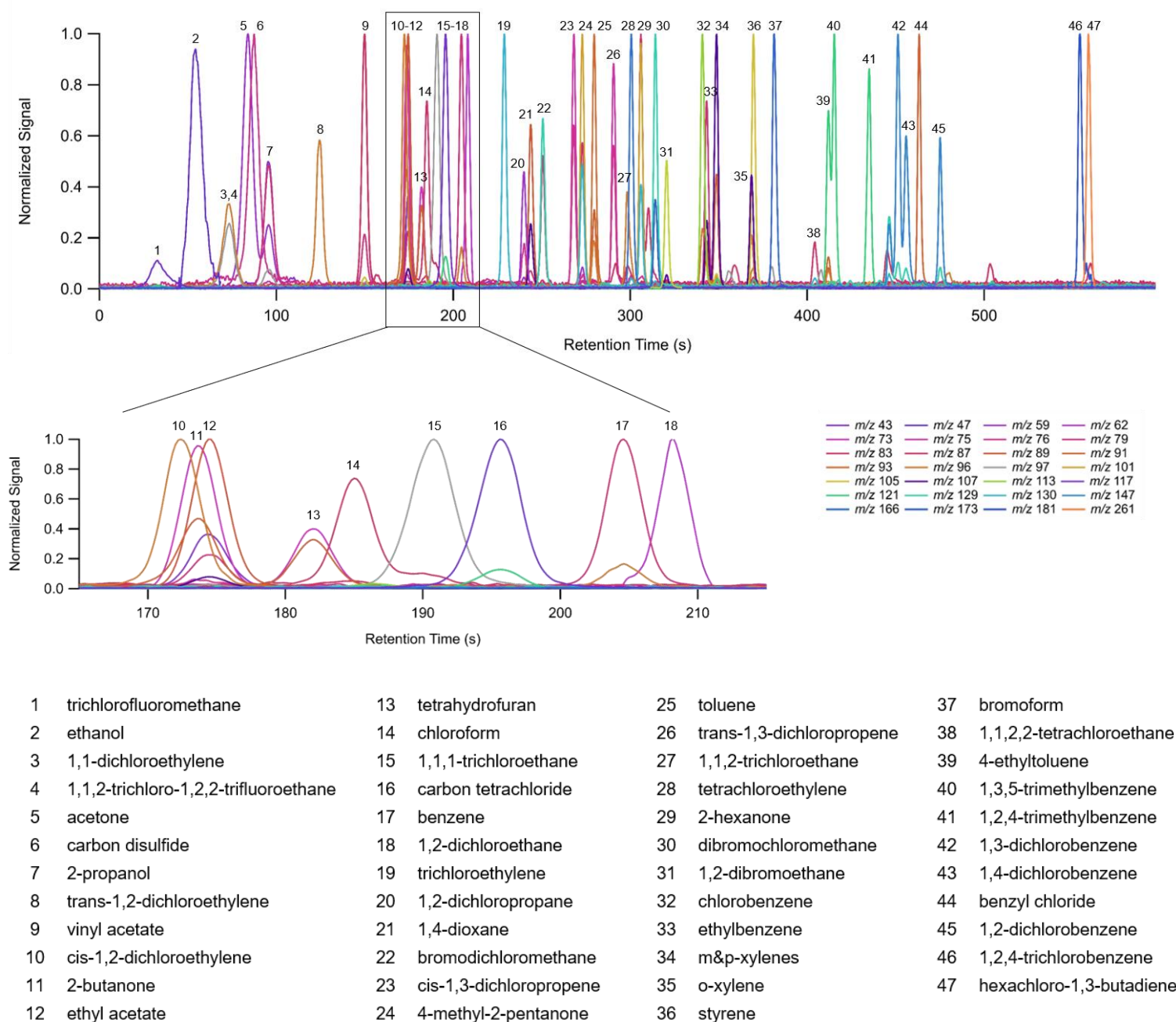




## Exemplary ARI TDPC-GC-EI-TOF Chromatogram

Aerodyne TDPC-GC chromatographic separation of a selection of atmospherically relevant VOCs, detected with EI-TOF-MS

Figure shows 10 minute chromatographic separation of 60 atmospheric VOCs. Note that species are separated into 4 plots for clarity, but are resolved in a single 10 minute chromatogram.

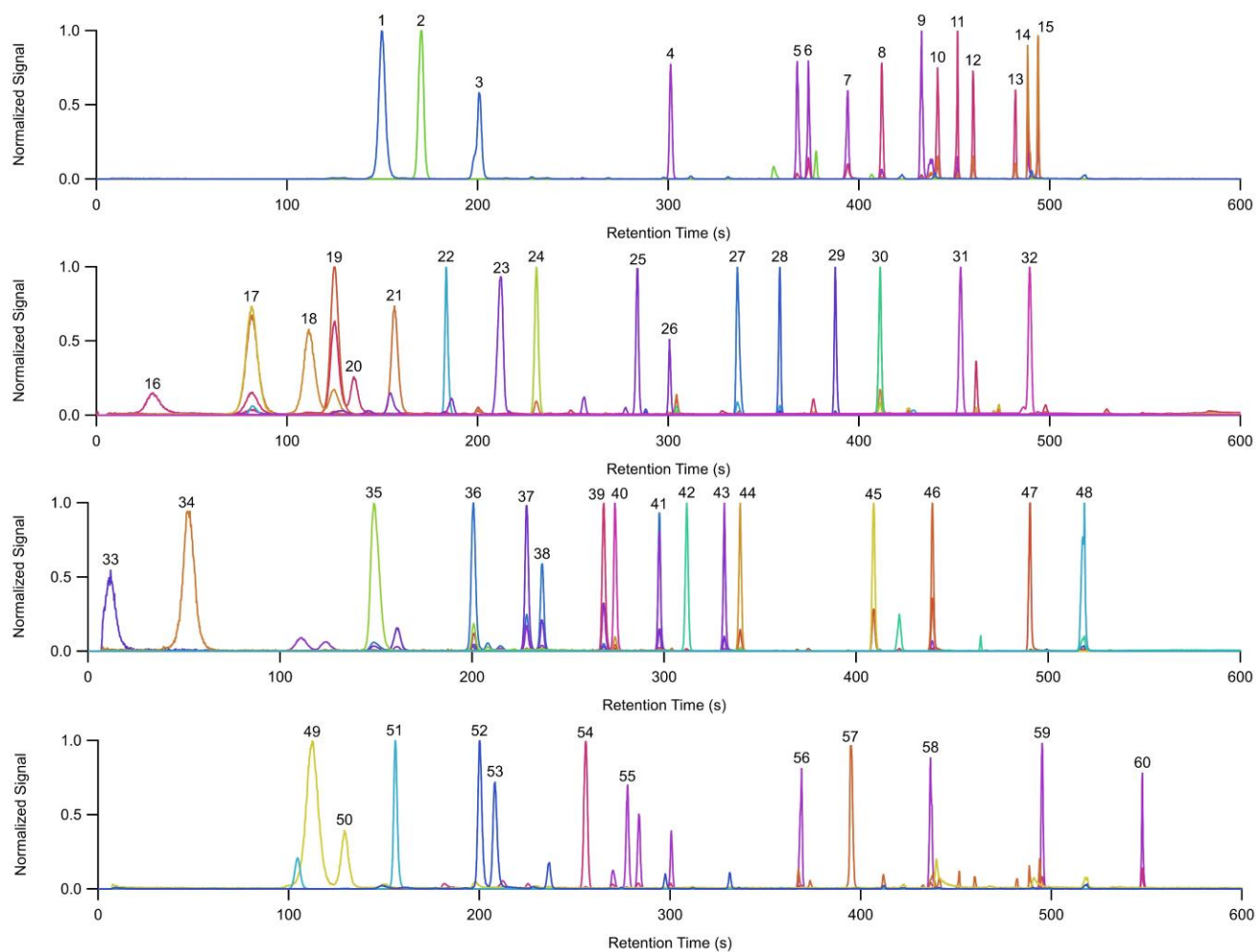

1 trans-C <sub>2</sub> H <sub>2</sub> Cl <sub>2</sub>	16 methanol	31 beta-pinene	46 4-ethyl toluene
2 CCl <sub>4</sub>	17 isoprene	32 cymene	47 benzyl chloride
3 ethyl acetate	18 acetone	33 butadiene	48 D5-siloxane
4 toluene	19 isopropanol	34 vinyl bromide	49 carbon disulfide
5 ethyl benzene	20 acetonitrile	35 methyl tert-butyl ether	50 allyl chloride
6 m&p-xylenes	21 acrylonitrile	36 vinyl acetate	51 n-hexane
7 o-xylene	22 tert-amyl ethyl ether	37 2,2,4-trimethyl pentane	52 methyl ethyl ketone
8 i-propyl benzene	23 cyclohexane	38 n-heptane	53 tetrahydrofuran
9 n-propyl benzene	24 benzene	39 1,4-dioxane	54 methyl cyclohexane
10 3-ethyl benzene	25 3-methyl heptane	40 CHCl <sub>2</sub> Br	55 2-methyl heptane
11 2-ethyl benzene	26 n-octane	41 4-methyl-2-pentanone	56 n-nonane
12 1,2,4-trimethyl benzene	27 D4-siloxane	42 D3-siloxane	57 styrene
13 1,3,5-trimethyl benzene	28 1,2,3-trimethyl benzene	43 2-hexanone	58 n-decane
14 m-diethyl benzene	29 limonene	44 CHBr <sub>2</sub> Cl	59 n-undecane
15 p-diethyl benzene	30 alpha-pinene	45 CHBr <sub>3</sub>	60 n-dodecane



## Linearities ( $R^2$ ) and limits of detection (LOD) for a variety of atmospherically relevant compounds measured by Aerodyne TDPC-GC-Vocus PTR-TOF-MS

These  $R^2$  and LOD values are measurements from 1 L pre-concentrated sample volumes collected at 100 sccm over analyte concentration ranges of 0.1 – 10 ppb.

Note that this table gives examples measured by the TDPC-GC-Vocus PTR-TOF-MS system, but compounds that can be resolved by the system **are not limited to** those specified on this list.

	$R^2$	LOD (ppt) <sup>a,b</sup>		$R^2$	LOD (ppt) <sup>a,b</sup>
<b>Aromatic Hydrocarbons</b>			<b>Halocarbons</b>		
1 benzene	0.998	1.5	21 bromoform	0.997	1.3
2 toluene	0.989	0.8	22 carbon tetrachloride	0.999	0.4
3 styrene	0.990	0.7	23 chloroform	1.000	0.7
4 ethyl benzene	0.997	1.0	24 bromodichloro methane	0.999	0.8
5 m&p-xylenes	0.979	0.6	25 dibromochloro methane	0.996	0.9
6 o-xylene	0.977	1.5	26 trichlorofluoro methane	0.996	4.1
7 4-ethyl toluene	0.965	2.2	27 1,1-dichloro ethylene	0.998	2.8
8 1,2,4-trimethyl benzene	0.997	1.3	28 cis-1,2-dichloro ethylene	1.000	2.5
9 1,3,5-trimethyl benzene	0.994	1.1	29 trans-1,2-dichloro ethylene	1.000	2.7
			30 tetrachloro ethylene	0.996	1.8
<b>Misc. Oxygenates</b>			31 trichloro ethylene	0.999	1.2
10 acetone	0.999	3.1	32 1,1,1-trichloro ethane	0.998	0.9
11 methyl ethyl ketone	0.993	1.4	33 1,1,2,2-tetrachloro ethane	0.997	2.2
12 2-hexanone	0.998	0.4	34 1,1,2-trichloro ethane	0.997	4.7
13 4-methyl-2-pentanone	0.997	0.3	35 1,2-dibromo ethane	0.998	0.8
14 ethanol	0.995	30	36 1,2-dichloro ethane	1.000	1.7
15 2-propanol	0.981	40	37 cis-1,3-dichloro propene	0.998	0.4
16 ethyl acetate	1.000	0.7	38 trans-1,3-dichloro propene	0.997	0.4
17 vinyl acetate	0.998	4.0	39 chloro benzene	0.997	0.8
18 1,4-dioxane	0.993	1.8	40 benzyl chloride	0.996	3.1
19 tetrahydro furan	0.998	2.5	41 1,2-dichloro benzene	0.996	7.2
20 carbon disulfide	0.999	1.6	42 1,3-dichloro benzene	0.992	1.6
			43 1,4-dichloro benzene	0.996	2.6

<sup>a</sup> Limit of detection reported in units of parts-per-trillion (ppt)

<sup>b</sup> Limit of detection calculated as  $3\sigma$  of the baseline, measured at  $3x$  FWHM before and after elution time.



## Linearities ( $R^2$ ) and limits of detection (LOD) for a variety of atmospherically relevant compounds measured by Aerodyne TDPC-GC-EI-TOF

These  $R^2$  and LOD values are measurements from 1 L pre-concentrated sample volumes collected at 100 sccm over analyte concentration ranges of 0.1 – 10 ppb.

Note that this table gives examples measured by the TDPC-GC-EI-TOF system, but compounds that can be resolved by the system **are not limited to** those specified on this list.

	$R^2$	LOD (ppt) <sup>a,b</sup>		$R^2$	LOD (ppt) <sup>a,b</sup>		
<b>Hydrocarbons</b>			<b>Halocarbons</b>				
1	isoprene	0.996	0.8	41	bromoform	0.996	0.1
2	cis-2-pentene	0.933	6.2	42	carbon tetrachloride	0.998	0.1
3	trans-2-pentene	0.957	6.0	43	chloroform	0.999	0.3
4	n-pentane	0.996	3.0	44	bromodichloro methane	0.999	0.2
5	i-pentane	0.980	5.0	45	dibromochloro methane	0.996	0.2
6	methyl cyclopentane	0.999	0.7	46	dichloro methane	1.000	0.9
7	1-hexene	1.000	6.6	47	trichlorofluoro methane	0.990	0.7
8	n-hexane	0.998	3.0	48	1,1-dichloro ethylene	1.000	0.3
9	cyclohexane	1.000	0.6	49	cis-1,2-dichloro ethylene	0.999	0.1
10	methyl cyclohexane	0.997	0.2	50	trans-1,2-dichloro ethylene	0.998	0.8
11	3-methyl pentane	0.998	2.4	51	tetrachloro ethylene	0.995	0.1
12	2-methyl hexane	0.993	0.3	52	trichloro ethylene	0.999	0.5
13	3-methyl hexane	0.979	0.7	53	1,1,1-trichloro ethane	0.998	0.2
14	n-heptane	0.998	0.8	54	1,1,2,2-tetrachloro ethane	0.997	0.6
15	2,3-dimethyl pentane	0.999	0.4	55	1,1,2-trichloro ethane	0.996	0.3
16	2,4-dimethyl pentane	0.999	2.1	56	1,1-dichloro ethane	1.000	0.3
17	2-methyl heptane	0.998	1.1	57	1,2-dibromo ethane	0.998	0.2
18	3-methyl heptane	0.999	1.0	58	1,2-dichloro ethane	1.000	0.5
19	n-octane	0.999	0.4	59	cis-1,3-dichloro propene	0.995	0.1
20	2,2,4-trimethyl pentane	0.999	0.1	60	trans-1,3-dichloro propene	0.995	0.3
21	2,3,4-trimethyl pentane	0.997	0.4	61	1,2-dichloro propane	0.999	0.3
22	alpha-pinene	0.991	0.2	62	chloro benzene	0.992	0.2
23	limonene	0.998	0.6	63	benzyl chloride	0.994	0.8
24	n-nonane	0.999	0.1	64	1,2-dichloro benzene	0.996	0.5
25	n-decane	1.000	0.2	65	1,3-dichloro benzene	0.998	0.2
26	n-undecane	1.000	0.5	66	1,4-dichloro benzene	0.995	0.4
27	n-dodecane	0.994	0.5				
				<b>Aromatic Hydrocarbons</b>			
				67	benzene	1.000	0.3
				68	toluene	0.996	0.2
				69	styrene	0.999	0.2
				70	ethyl benzene	0.998	0.2
				71	m&p-xylenes	0.999	0.3
				72	o-xylene	0.997	0.5
				73	i-propyl benzene	0.997	0.3
				74	n-propyl benzene	0.997	0.3
				75	2-ethyl toluene	0.997	1.3
				76	3-ethyl toluene	0.996	0.6
				77	4-ethyl toluene	0.996	0.6
				78	1,2,3-trimethyl benzene	0.996	0.2
				79	1,2,4-trimethyl benzene	0.997	0.6
				80	1,3,5-trimethyl benzene	0.999	0.5
				81	1,3-diethyl benzene	0.997	0.3
				82	1,4-diethyl benzene	0.998	0.3
<b>Misc. Oxygenates</b>							
28	acetone	1.000	1.0				
29	methyl ethyl ketone	0.999	0.9				
30	2-hexanone	1.000	0.5				
31	4-methyl-2-pentanone	1.000	0.6				
32	ethanol	1.000	1.2				
33	2-propanol	1.000	0.2				
34	ethyl acetate	0.996	0.8				
35	vinyl acetate	0.999	0.5				
36	1,4-dioxane	0.999	0.3				
37	methyl tert-butyl ether	1.000	0.8				
38	amyl ethyl ether	0.999	0.2				
39	tetrahydro furan	0.999	0.9				
40	carbon disulfide	1.000	0.7				

<sup>a</sup> Limit of detection reported in units of parts-per-trillion (ppt)

<sup>b</sup> Limit of detection calculated as  $3\sigma$  of the baseline, measured at 3x FWHM before and after elution time.