

OpenView Series

The universal camera for laser beam analysis and imaging

The largest spectral range

The **OpenView** is an universal camera offering the largest spectral range from UV to THz domain. Our technology is based on a high performance photon to IR converter able to convert any photon from 0,1 to 3 000 µm.

Users select it's own spectral band of the OpenView using optical filters in UV, Visible, IR, THz region. And take also benefits of a large sensitivity area ($\phi = 50$ mm) for high-power laser profiling and imaging.

Specifications		
Spectral range	Select your own spectral range	
Number of pixel	320x256	640x512
Active area (mm)	54,4 x 43,5	51,2 x 41
Spatial resolution	170 µm	80 µm
Minimum signal detection	50 μW/cm²	
Damage threshold	1 W / cm²	
Included Software	Vision and acquisition	
Prodcuts size (mm)	90 x 90 x 200	
Working temperature	Room temperature	
Supply Voltage	110/220 V	
Plug-in	Gigabit Ethernet	

Applications

Laser beam analysis (profiling, M2, divergency, ...)

2D or 3D multispectral imaging

Non destructive testing and industrial vision (composit, wood, plastic, ceramics, ...)

OpenView Series



OpenViewUV

Ö,05 - 1 μm









OpenViewIR 1,1 - 25 µm


OpenViewTHz 25 - 3000 µm

OpenViewMS 0,1 - 3000 µm


Vison and acquisition software included




UV Laser $\lambda = 0.2 \mu m$

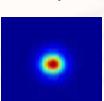

IR Laser THz Laser $\lambda = 3 \mu m$ $\lambda = 2800 \, \mu m$

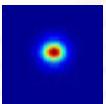

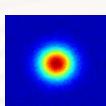

Key benefits

Choice of spectral range: UV, Visible, IR, THz or multispectral (0,1 to 3000 µm)

Largest and uncooled detection surface

Adapted for all high-power lasers sources (Excimer, solid state, OPA, QCL, CO2, Gunn diode, BWO, .)

Noise reducing software

Accessories



OPFilter-VIS



OPFilter-THz



OPFilter-IR



OPFilter-UV






OpenView Series

The universal camera for laser bean analysis and imaging


View and characterize your beam

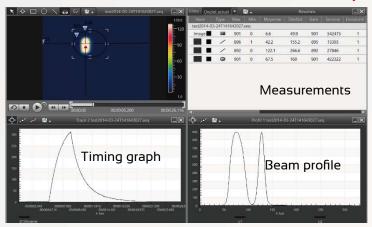

Functionnality

Real time vision/acquisition mod Profile analysis Timing graph and measurements

Noise and Thermal software processing to enable low-energy

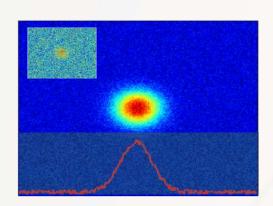


OPA Laser

Spectral range = 0,1 - 3 THz
Peak energy : 1 µJ
Beam diameter : 0,9 mm
Pulse duration : 100 fs
Frenquency rate : 1 kHz

Increase by 10 to 100 times the signal to noise ratio

Unique THz beam profiling characterization

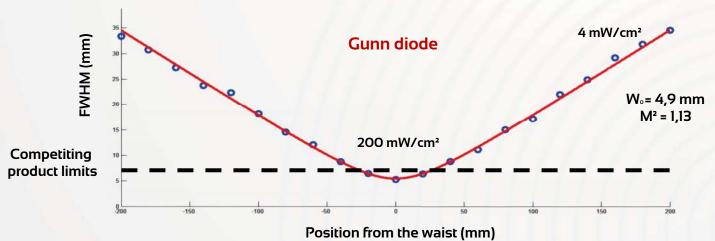



NETHIS - NEW TERAHERTZ IMAGING SYSTEMS

Tel. + (33) 5 47 74 62 10

More information: info@nethis-thz.com