

Thermo Scientific ARL EQUINOX 5000 Series

High Resolution X-Ray Diffractometer

The Thermo Scientific™ ARL™ EQUINOX 5000 is a multi purpose X-ray diffractometer designed specifically for high level research applications including powder analysis, thin layer analysis, texture, stress, and micro diffraction.

The system utilizes a unique tall vertical goniometer in θ/θ or $\theta/2\theta$ mode and the user has a large choice of features including higher resolution detectors and a large choice of X-ray sources. The free space in the sample compartment is very large and will accommodate practically any sample and any sample handling device including furnaces, auto samplers, and large assemblies.



Perfect for advanced research applications.



Research

Environment / Geology



Polymers /
Biomaterials


Metallurgy



Electronics

Photovoltaic

The EQUINOX 5500 is an evolution of EQUINOX 5000 model. It is exactly the same configuration except that it works with a bigger radius of curvature detector. So, obtained data have a better resolution.

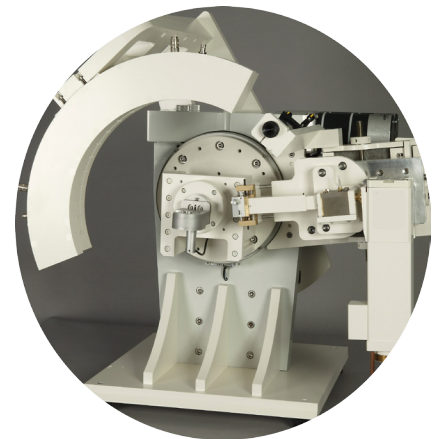


Real time acquisition
Extremely accurate
Easy to use
No alignment needed
Very reliable
Superb resolution
No maintenance


Thermo
SCIENTIFIC

Specifications

X-ray source	Generator: 3500 W (60kV/60mA - Option 30kV limitation Standard sealed X-Ray tube (Cu, Co, Cr, Mo, ... radiation)
Detector	Curved Position Sensitive X-ray Detector, CPS120 Curvature radius: 250mm Acquisition in real time over 120° 2Theta Option: CPS590 (500mm, 90° 2Theta – EQUINOX 5500)
Goniometric table	Theta / Theta or Theta / 2 theta column Motorization per step with encoder, up to 5/10000° Acquisition in asymmetric mode
Sample holders	Fixed, with or without sample rotation Reflection or transmission mode Sample changer: 30 positions Thin film attachment (reflectometry, grazing incidence) Chi/Phi attachment (texture, stress) Mapping attachment (X, Y, Z translations)
In-situ camera	Heating camera: low and high temperature Humidity chamber, Electrochemistry cell Tensile stage
Optic	Monochromatic
Computer	Windows© 8, 7 Vista, XP
Software	Real time diffractogram display Multiple and automatic recording Peak search and deconvolution with several shapes Search match and quantification Degree of crystallinity determination Cell parameters, crystallite size, lattice strain Crystal structure analysis Microstructure analysis Rietveld analysis Phase transition, dynamic studies Thin film and multilayer analysis Free open database, COD Option: ICDD PDF2 or PDF4 databases
Power	30-32A / 208-230V / 50-60 Hz
Cooling water	Flow: ≈ 3,5 L/min - T°C: ≈18 to 25 °C
Dimensions	1815 mm H x 744 mm D x 1206 mm W EQUINOX 5000 1966 mm H x 1445 mm D x 1433 mm W EQUINOX 5500

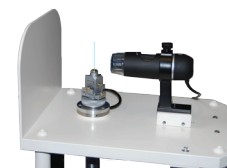

You want to perform transmission analysis? Choose either the Sample Spinner for reflection and transmission analysis on Powder samples or the Capillary sample holder for Transmission mode analysis depending on your needs.

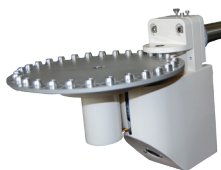

You need to analysis several powder samples a day, to reduce dramatically operator's time? Optimize your throughput with the 30-positions sample changer.




You need to perform grazing incidence measurement? Try our Thin Layer analysis attachment.



HTK1200N

This system allows you to work with controlled temperature or atmosphere chamber. Anton Paar and other suppliers chambers can be available depending on your needs.

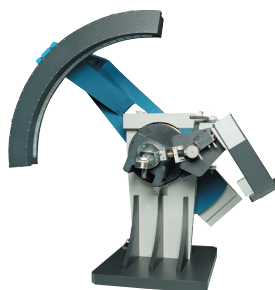

High resolution diffractometer with a CPS 590 curved detector.

www.thermoscientific.com/xrd

©2016 Thermo Fisher Scientific Inc. All rights reserved. Microsoft and Windows are registered trademarks of Microsoft Corporation. ISO is a trademark of the International Standards Organization. All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.



Thermo Fisher Scientific (Ecuiblen) SARL, Switzerland is ISO certified.



Africa +43 1 333 50 34 0
Australia +61 3 9757 4300
Austria +43 810 282 206
Belgium +32 53 73 42 41
Canada +1 800 530 8447
China +86 21 6865 4588

Denmark +45 70 23 62 60
Europe-Other +43 1 333 50 34 0
Finland/Norway/Sweden
+46 8 556 468 00
France +33 2 38 80 45 45
Germany +49 6103 408 1014

India +91 22 6742 9494
Italy +39 02 950 591
Japan +81 45 453 9100
Latin America +1 561 688 8700
Middle East +43 1 333 50 34 0
Netherlands +31 76 579 55 55

New Zealand +64 9 980 6700
Russia/CIS +43 1 333 50 34 0
Spain +34 914 845 965
Switzerland +41 61 716 77 00
UK +44 1442 233555
USA +1 800 532 4752

Thermo
SCIENTIFIC

A Thermo Fisher Scientific Brand